

Instructions for use Part 1

Oticon Engage
miniRITE



Introduction to this booklet

This booklet gives you guidance on how to use and maintain your new hearing aid. Please read the booklet carefully, including the Warnings section. This will help you to get the most out of your new hearing aid.

Your hearing care professional has adjusted the hearing aid to meet your needs. If you have additional questions, contact your hearing care professional.


A hearing care professional* (hearing aid professional, audiologist, ENT (ear, nose and throat) doctor, and hearing aid dispenser) is a person who is appropriately educated and has proven competency in professionally assessing hearing, selecting, fitting, and delivering hearing instruments and rehabilitation care to persons with hearing loss.

The education of the hearing care professional is in accordance with national or regional regulations.


FW 6.1

Oticon Engage miniRITE
GTIN: (01) 05714464013234

Your hearing aid, speaker, and earpiece miniRITE

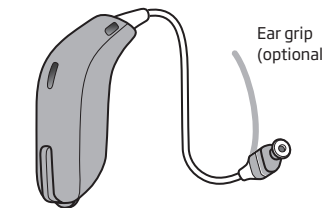


The hearing aid uses one of the following speakers:

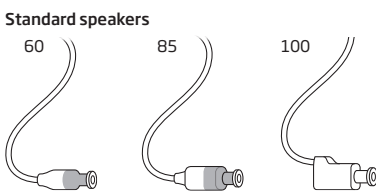

Power flex mould speakers


- ☐ 100 (Power Instrument)
- ☐ 105 (Power Instrument)

Dome sizes

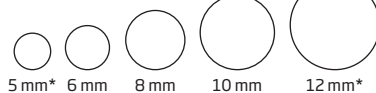

* only as open dome for speaker 60
** not for open dome

Hearing aid What it is and does

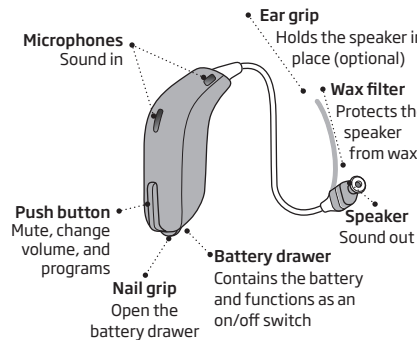

The speaker uses one of the following earpieces:

Standard earpieces

- ☐ Open dome
- ☐ Bass dome, single vent
- ☐ Bass dome, double vent
- ☐ Power dome
- ☐ Grip Tip

Available in small and large, left and right, with or without vent.


Customised earpieces

- ☐ Micro mould
 - ☐ LiteTip
 - ☐ VarioTherm® micro mould
 - ☐ VarioTherm® LiteTip
- ® VarioTherm is a registered trademark of Dreve


Turn hearing aid ON and OFF

The battery drawer is used to switch the hearing aid on and off. To save battery life, make sure your hearing aid is switched off when you are not wearing it. If you wish to return to the standard settings of the hearing aid, simply open and close the battery drawer (quick reset).

Turn ON
Close the battery drawer with the battery in place.

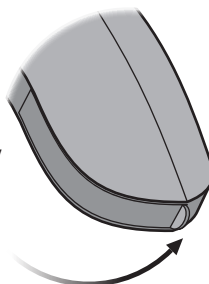

Turn OFF
Open the battery drawer.

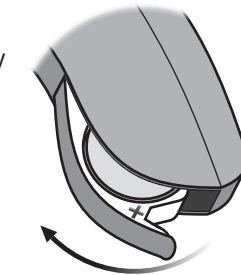

When to replace the battery

When it is time to replace the battery, you will hear three tones repeated in moderate intervals until the battery runs out.

Three alternate tones*
= The battery is running low

Four descending tones
= The battery has run out

Battery tip

To make sure the hearing aid is always working, bring spare batteries with you, or replace the battery before you leave home.

Optional LED


Continuous red flashes indicate low battery.

Note: Batteries need to be replaced more often if you are streaming to your hearing aid.

* Bluetooth® will be turned off and it will not be possible to use wireless accessories.

How to replace the battery (size 312)

1. Remove
Fully open the battery drawer.
Remove the battery.

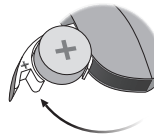


2. Uncover
Remove the sticky label from the + side of the new battery.




Tip:
Wait 2 minutes so that the battery can draw air, to ensure optimal functioning.

3. Insert
Insert the new battery into the battery drawer with the + side facing upwards.

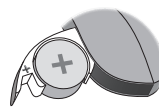


4. Close
Close the battery drawer.
The hearing aid will play a jingle through the earpiece.

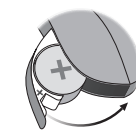

Hold the earpiece close to your ear to hear the jingle.

Tip
The MultiTool can be used to change the battery. Use the magnetic end to remove and insert batteries.

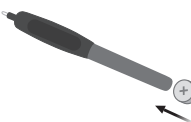


The MultiTool is provided by your hearing care professional.

Identify left and right hearing aid

It is important to distinguish between the left and the right hearing aid, as they may be programmed differently. You can find left/right colour indicators on the hearing aid itself and on 60 and 85 speakers as shown. Indicators (either L or R) can also be found on 100 speakers and some earpieces.

A RED indicator marks the RIGHT hearing aid.

A BLUE indicator marks the LEFT hearing aid.

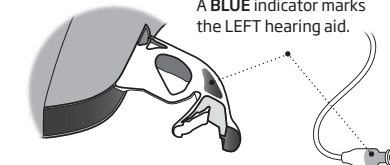

Put on hearing aid


Step 1

Place the hearing aid behind your ear. The speaker should always be used with an earpiece attached. Use only parts designed for your hearing aid.



Step 2

Hold the bend of the speaker wire between your thumb and index finger. The earpiece should point toward the opening of the ear canal.

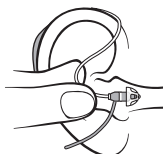


Step 3

Gently push the earpiece into your ear canal until the speaker wire sits close to your head. If the speaker has an ear grip, place it in the ear so it follows the contour of the ear.



Flight mode

When flight mode is activated, Bluetooth® is turned off. The hearing aid will still be working. Pressing the push button on one hearing aid will activate flight mode on both hearing aids.


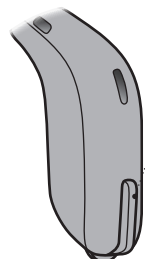


To activate and deactivate
Press the push button for at least 7 seconds. A jingle confirms your action.

Opening and closing the battery drawer will also deactivate flight mode.

Change volume

The push button allows you to adjust the volume. You will hear a click when you increase or decrease the volume.

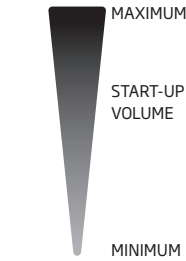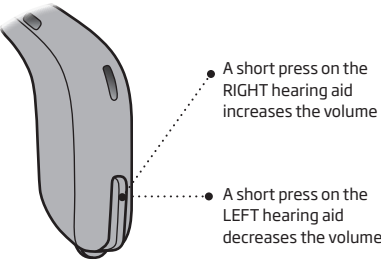


To be filled out by the hearing care professional

Volume change	☐ LEFT	☐ RIGHT	☐ Short press
---------------	-------------------------------	--------------------------------	--------------------------------------

Change program


Your hearing aid can have up to four different programs. These are programmed by your hearing care professional. You will hear one to four tones when you change program depending on the program.

See the Sound and LED indicator section.

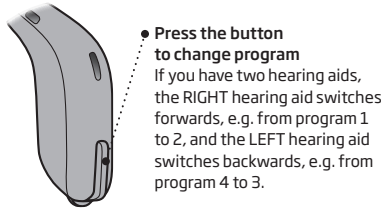

Mute

Use the mute function if you need to silence the hearing aid. The mute function only mutes the microphone(s) on the hearing aid.

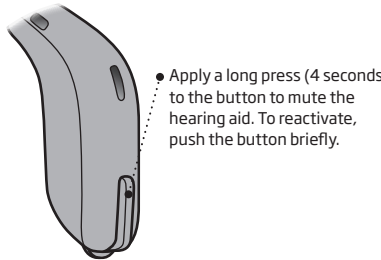

Use hearing aid with iPhone and iPad

Your hearing aid is Made for iPhone® and allows for direct communication and control with iPhone, iPad® or iPod touch®. For assistance in using these products with your hearing aid, please contact your hearing care professional.



To see how to pair you hearing aid with your iPhone please visit www.oticon.global/pairing

For information on compatibility, please visit www.oticon.global/compatibility

Use of the Made for Apple badge means that an accessory has been designed to connect specifically to the Apple product(s) identified in the badge, and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Please note that use of this Accessory with iPod, iPhone, or iPad may affect wireless performance.

Wireless accessories and other options

As an enhancement to your wireless hearing aid, a range of accessories is available. These can enable you to hear and communicate better in many everyday situations.

ConnectClip

When ConnectClip is paired with your mobile phone, you can use the hearing aid as a hands-free headset. ConnectClip can also be used as a remote microphone

TV Adapter 3.0

TV Adapter streams sound directly from a TV or electronic audio device to your hearing aid.

Remote Control 3.0

Offers the ability to change program, adjust volume, or mute your hearing aid.

Oticon ON app

An application that lets you control your hearing aid from your mobile phone or tablet. For iPhone, iPad, iPod touch, and Android devices. Ensure that you only download and install Oticon ON app from the official app stores.

Telecoil program

Telecoil helps you hear better when using a telephone with a built-in loop or when you are in buildings with teleloop systems such as theatres, churches, or lecture rooms. This symbol or a similar sign is shown wherever a teleloop has been installed.



Autophone - optional for Oticon Engage and Oticon Siya

Activates a phone program in the hearing aid, if you place a dedicated magnet on your phone.



CROS - optional for Oticon Engage, Oticon Opn Play 1, Oticon Ruby 1 and Ruby 2

Solution for people with an unaidable hearing loss in one ear. Oticon CROS on the poorer ear transmits sound to the hearing aid on the better ear.

For more information visit www.oticon.global or contact your hearing care professional.

Sound and LED indicators

Different sounds and LED lights indicate the hearing aid status.

The different indicators are listed on the following pages. For charger LED indicators, see Charger Instruction for use.

Your hearing care professional can set sound and LED indicators to match your preferences.

Program	☐ Sound	When to use
1	1 tone	
2	2 tones	
3	3 tones	
4	4 tones	

ON/OFF	Sound
On	☐ Jingle
Off	☐ 4 descending tones
Volume	Sound
Start-up volume	☐ 2 beeps
Minimum/maximum volume	☐ 3 beeps
Volume up/down	☐ 1 beep
Mute activated on ON app, ConnectClip or Remote Control 3.0	

Accessories	☐ Sound
TV Adapter	2 different tones
ConnectClip remote microphone	2 different tones
Flight mode	Sound
Flight mode active	4 descending tones + short jingle
Flight mode inactive	4 descending tones + short jingle

Warnings	Sound
Low battery	☐ 3 alternate tones
Battery shut down	4 descending tones
Microphone service check needed	8 beeps repated 4 times
No light in the hearing aid LED when placed in the charger	
The hearing aid LED flashes red/orange when the hearing aid is placed in the charger	

Your individual hearing aid settings

To be filled out by your hearing care professional.

Tinnitus SoundSupport: Limitation on use			
☐	No limitation on use		
☐	Program	Start-up volume (Tinnitus)	Max volume (Tinnitus)
☐	1	Max _____ hours per day	Max _____ hours per day
☐	2	Max _____ hours per day	Max _____ hours per day
☐	3	Max _____ hours per day	Max _____ hours per day
☐	4	Max _____ hours per day	Max _____ hours per day

Settings overview for your hearing aid			
Left		Right	
☐ Yes	☐ No	Change volume	☐ Yes ☐ No
☐ Yes	☐ No	Change program	☐ Yes ☐ No
☐ Yes	☐ No	Tinnitus SoundSupport	☐ Yes ☐ No
Volume control indicators			
☐ On	☐ Off	Beeps at min/max volume	☐ On ☐ Off
☐ On	☐ Off	Beeps when changing volume	☐ On ☐ Off
☐ On	☐ Off	Beeps at start-up volume	☐ On ☐ Off
Battery indicators			
☐ On	☐ Off	Low battery warning	☐ On ☐ Off

IMPORTANT NOTICE

Do not use the mute function as an off switch, as the hearing aids are still using battery power in this mode.