

Instructions for use Part 1

Oticon Sage
miniRITE T



Made for
iPhone | iPad | iPod

Works with
android

oticon

Model overview

This booklet is valid for the Oticon Sage family for the following hearing aid model:

FW 1.0

Oticon Sage miniRITE T
GTIN: (01) 05714464127238

Introduction to this booklet


This booklet guides you on how to use and maintain your new hearing aids. Ensure you read this booklet carefully, including the Warnings section. This will help you get the most benefit from your hearing aids. Your hearing care professional has adjusted the hearing aid to meet your needs. If you have

additional questions, contact your hearing care professional.


A hearing care professional* (hearing aid professional, audiologist, ENT (ear, nose and throat) doctor, and hearing aid dispenser) is a person who is appropriately educated and has proven competency in professionally assessing hearing, selecting, fitting, and delivering hearing instruments and rehabilitation care to persons with hearing loss. The education of the hearing care professional is in accordance with national or regional regulations.

* The job title may vary from country to country.

Your hearing aid, speaker, and earpiece

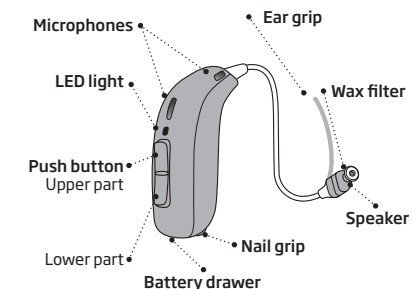

The hearing aid uses one of the following speakers:

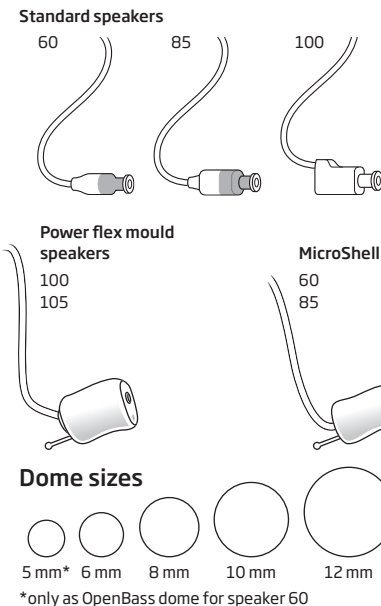

The speaker uses one of the following earpieces:

Standard earpieces

- ☐ OpenBass dome
 - ☐ Bass dome, double vent
 - ☐ Power dome
 - ☐ Grip Tip
- Available in small and large, left and right, with or without vent.

Customised earpieces

- ☐ Micro mould
 - ☐ LiteTip
 - ☐ VarioTherm® Micro mould
 - ☐ VarioTherm® LiteTip
- ® VarioTherm is a registered trademark of Dreve

Note

For details on replacing the dome, see the Replace standard earpieces section in instruction for use Part 2.

Battery performance

The battery performance varies depending on your individual use and hearing aid settings. Streaming sound from a TV, mobile phone or connectivity devices can influence this performance.


Turn hearing aids ON/OFF

Using the battery drawer

The battery drawer is used to turn the hearing aids ON and OFF. To save battery life, make sure your hearing aids are switched OFF when you are not wearing them. To perform a quick reset of hearing aid settings, open and close the battery drawer.


Turn ON

Close the battery drawer with the battery in place.

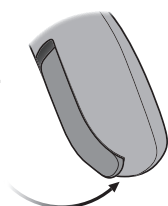

Turn OFF

Open the battery drawer.

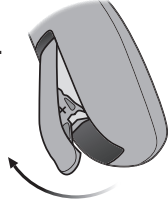

Identify left and right hearing aid

It is important to distinguish between the left and the right hearing aid, as they may be programmed differently. You can find left/right colour indicators on the hearing aid itself and on 60 and 85 speakers as shown.

Indicator markings (either L or R) can also be found on 100 speakers and customised earpieces.

For 105 speakers, the indicator is found on the earpiece.

Low battery indication

Just before the battery runs out completely, you hear four descending tones. To extend battery performance, ensure you stop any audio streaming.

When the battery is running low, you hear three alternate beeps. This gives you approximately 15 minutes before the hearing aid runs out of battery. At this point, Bluetooth® connectivity is turned OFF.

Three alternate beeps


Three beeps = The battery is running low

Four descending beeps

Four beeps = The battery has run out

LED light


Continuous ORANGE blinks indicate low battery.

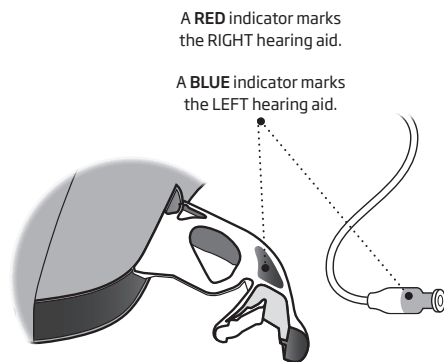

How to replace the disposable battery – size 312


1. Remove

Fully open the battery drawer and remove the battery.

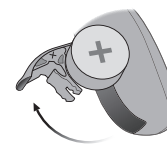


2. Uncover

Remove the sticky label from the + side of the new battery.

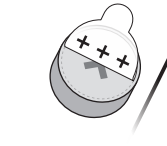


3. Insert

Insert the new battery into the battery drawer with the + side facing upwards.

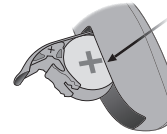

4. Close

Close the battery drawer. You may hear a jingle through the earpiece.

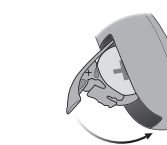

Put on hearing aid

Step 1


Place the hearing aid behind your ear. You should always use the speaker with an earpiece attached.



Ensure you only use parts designed for your hearing aid.


Step 2

Hold the bend of the speaker wire between your thumb and index finger. The earpiece should point towards the opening of the ear canal.

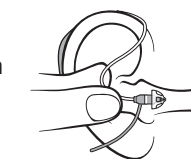


Step 3

Gently push the earpiece into your ear canal until the speaker wire sits close to your head. If the speaker has an ear grip, place it in the ear so it follows the contour of the ear.

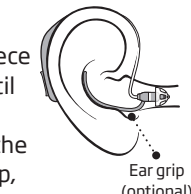

Flight mode

When Flight mode is activated, Bluetooth connectivity is turned OFF. However, the hearing aid is still turned ON and functioning. Be aware that pressing the push-button on one hearing aid, activates Flight mode on both hearing aids. For more information about sounds and lights, see the Sound and LED light indicators section.

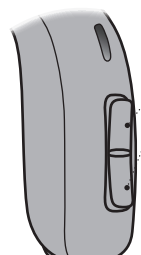


To activate and deactivate
Press and hold either end of the push-button for at least seven seconds. A jingle and an LED light pattern confirm your action.

Opening and closing the battery drawer also deactivates Flight mode.


Change volume

The push-button lets you adjust the volume. When you increase or decrease the volume, you hear a beep.

For information regarding button press times, see the table General settings overview for your hearing aid, in the Your individual hearing aid settings section at the end of this booklet.

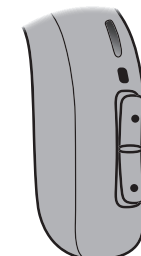

- Press button to increase volume
- Press button to decrease volume

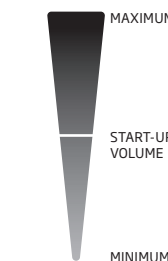

000026683000001


oticon

Change program

Your hearing aid can save up to four different programs configured by your hearing care professional. Depending on the program you choose (1,2,3 or 4), you hear one to four tones.

For information regarding tones, see the Sound and LED light indicators section.

For information regarding button press times see the table General settings overview for your hearing aid, in the Your individual hearing aid settings section at the end of this booklet.

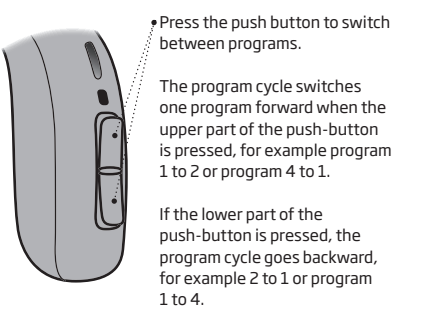

Mute your hearing aids

You can mute your hearing aids by using one of the following optional devices/app:

- Oticon Companion
- ConnectClip
- Remote Control 3.0

How to unmute your hearing aids

You can unmute your hearing aids by using one of the optional devices/app or

by applying a short press to the upper or lower part of the push-button on the hearing aids.

You can also mute your hearing aids by pressing either end of the push-button for four seconds.

Use your hearing aid with iPhone, iPad and iPod

Your hearing aids are Made for iPhone hearing aids and allow for direct streaming from your iPhone, iPad or iPod touch.* Oticon Companion can be used to control your hearing aids from your mobile device.* Additionally, your hearing aids offer hands-free communication with compatible iPhone and iPad devices.*

For assistance in using your hearing aids with any of these products, contact your hearing care professional.

Use your hearing aid with Android devices

Your hearing aids support Audio Streaming for Hearing Aids (ASHA) and allow for direct streaming from selected Android™ devices.* Oticon Companion can be used to control your hearing aids from your mobile device.*

Apple, the Apple logo, iPhone, iPad, and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc. Android, Google Play, and the Google Play logo are trademarks of Google LLC. Use of the Made for Apple badge means that an accessory has been designed to connect specifically to the Apple product(s) identified in the badge, and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Android, Google Play, and the Google Play logo are trademarks of Google LLC.

For assistance in using your hearing aids with Android devices, contact your hearing care professional.


Pairing and compatibility

For instructions on how to pair your hearing aid with iPhone, iPad, iPod touch or Android devices, visit: www.oticon.global/pairing

*For a list of compatible iPhone, iPad, iPod touch and Android devices, visit: www.oticon.global/compatibility

Call handling

You can answer, reject or end phone calls with the push-buttons on your hearing aids. For this functionality, your hearing aids must be paired with a compatible iPhone*.

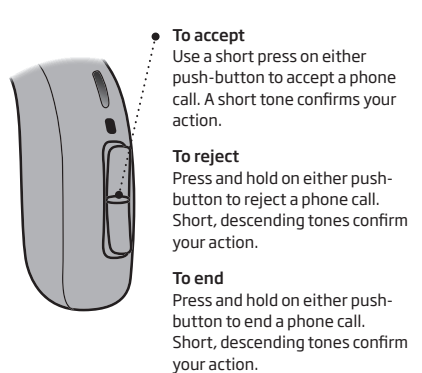

* For a list of compatible iPhones, visit: www.oticon.global/compatibility

Wireless accessories and other options

There are a range of accessories available as an enhancement to your wireless hearing aid. These enable you to hear and communicate better in everyday situations.

ConnectClip

A device that can be used as remote microphone and hands-free headset when paired to your mobile phone.

TV Adapter 3.0

A device that streams sound from a TV or electronic audio device, to your hearing aids.

Remote Control 3.0

A device that lets you change program, adjust volume, or mute your hearing aids.

Oticon Companion

An application that gives you wireless access to your hearing aids through a mobile device for control of functionality, and that allows an update of your hearing aids remotely. For iPhone, iPad, iPod touch and Android devices. Ensure that you only download Oticon Companion from the official app stores.

Phone Adapter 2.0

A device that when used together with hearing aids and ConnectClip, lets you communicate hands-free via a landline phone.

 **Telecoil**
Telecoil can help you hear better when using a phone with a built-in loop or when in buildings with teleloop systems such as theatres, places of worship, or lecture rooms. This symbol is shown wherever a teleloop has been installed.

EduMic

A device that can be used as a remote microphone in classrooms, work environments, public places (using Telecoil), and other settings.

Sound and LED light indicators

Different sounds and LED lights indicate the hearing aid status.

The different indicators are listed on the following pages.

Program	☐ Sound	☐ LED light*	When to use
1	1 tone		
2	2 tones		
3	3 tones		
4	4 tones		

*LED light blinks continuously or is repeated three times with short pauses

ON/OFF	Sound	LED light	LED light comments
ON	☐ jingle		Shown once
Volume	Sound	LED light	
Start-up volume	☐ 2 beeps		
Minimum/maximum volume	☐ 3 beeps		
Volume up/down	☐ 1 beep		
Mute			
			Continuous or repeated three times

Accessories	☐ Sound	☐ LED light	LED light comments
TV Adapter 3.0	2 different tones		Continuous or repeated three times
ConnectClip	2 different tones		
Flight mode	☐ Sound	☐ LED light	
Flight mode activated	Short jingle		
Flight mode deactivated	Short jingle	*	

Warnings	Sound	LED light	LED light comments
Low battery	☐ 3 alternate beeps		Continuously blinking
Battery shut down	4 descending tones		
Microphone service check needed	8 beeps repeated 4 times		Repeated four times

Long GREEN blink Short GREEN blink Long ORANGE blink Short ORANGE blink

IMPORTANT NOTICE
Do not use the mute function as an OFF button, as the hearing aid is still using battery power in this mode.

Your individual hearing aid settings

To be filled out by your hearing care professional.

General settings overview for your hearing aid					
Left			Right		
☐ Yes	☐ No	Change volume	☐ Yes	☐ No	
☐ Yes	☐ No		☐ Yes	☐ No	
☐ Short press		Change program	☐ Short press		
☐ Long press			☐ Long press		
Volume control indicators					
☐ ON	☐ OFF	Beeps at min/max volume	☐ ON	☐ OFF	
☐ ON	☐ OFF	Beeps when changing volume	☐ ON	☐ OFF	
☐ ON	☐ OFF	Beeps at start-up volume	☐ ON	☐ OFF	
Battery indicators					
☐ ON	☐ OFF	Low battery warning	☐ ON	☐ OFF	

To be filled out by your hearing care professional.

Tinnitus SoundSupport: Limitation on use			
☐	No limitation on use		
	Program	Start-up volume (Tinnitus)	Max volume (Tinnitus)
☐	1	Max _____ hours per day	Max _____ hours per day
☐	2	Max _____ hours per day	Max _____ hours per day
☐	3	Max _____ hours per day	Max _____ hours per day
☐	4	Max _____ hours per day	Max _____ hours per day
Tinnitus SoundSupport settings overview for your hearing aid			
Left		Right	
☐ Yes	☐ No	Tinnitus SoundSupport	☐ Yes ☐ No
☐ A) How to change Tinnitus SoundSupport volume in each ear separately To increase or decrease the volume (on one hearing aid only), use a short press on the upper or lower part of the push-button repeatedly until you reach your desired level.			
☐ B) How to change Tinnitus SoundSupport volume in both ears simultaneously You can use one hearing aid to increase/decrease the sound in both hearing aids. When changing the volume in one hearing aid, the volume on the other hearing aid follows. To increase volume, use a short press on the upper part of the push-button repeatedly. To decrease volume, use a short press on the lower part of the push-button repeatedly.			