


6 reasons why ReSound Smart Fit™ will enrich your fitting experience

ReSound GN



Download now on gnhearing.com

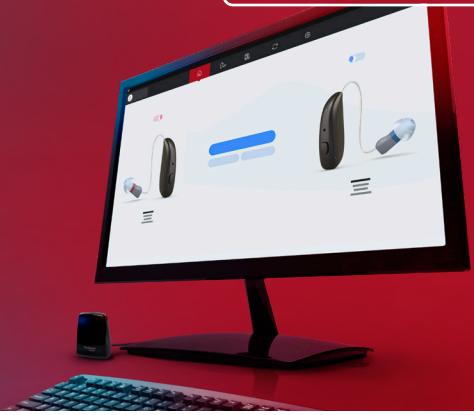

As a hearing care professional, we know you spend a lot of your day in fitting software. And you deserve the best experience. So, we're enriching your fitting experience, to ensure you can get more intuitive fittings in even fewer clicks.



Show – don't just tell

Introducing our new range of RIE (Receiver-in-Ear) demo styles that lets you change technology levels easily with 2 clicks¹, so you can let your clients hear and experience the difference advanced sound features make.



Spend less time clicking

We've made a 33-second improvement to click-through speed² and deliver a 62% click reduction¹ when clicking around for preferences and settings.



Personalised onboarding for your clients

Clients getting overwhelmed with information? Help your clients by printing out personalised onboarding reports with all the features and guidance they need.



Enjoy a more intuitive experience

With a significant simplification of user interface, there's now 74% less clutter² on the home screen. Your journey through our software has never been cleaner and simpler.



Revamped home screen

You can access software updates, counseling tools and learning resources from the home screen, making sure you're always up to date on software, accessory firmware, and all our latest materials.



Compare our different hearing aids and technology levels

Now, you can not only look up and check all our hearing aid models, but you can also compare the different devices across technology levels.

We've prioritised the fundamentals

Over many years, we've kept improving our fitting software. Every advanced feature you experience is built on a solid foundation of consistency and we are not stopping with the latest version of ReSound Smart Fit. Like any fine-tuning, we'll keep making it even better.

**"I think it is visually very appealing.
It's larger, it's easier to read, it's less
cumbersome with all the stuff popping up"**

Hearing Aid Fitter, 7 years experience, USA

**"Sometimes you need a quick reference
on what the serial number is. It's nice
that's available here"**





Hearing Aid Fitter, 1 years experience, USA

¹ GN Proprietary data on file

² Resound Smart Fit 2.0 compared to ReSound Smart Fit 1.17, technical data on file

6 easy steps to fit ReSound hearing aids

- 1 Launch ReSound Smart Fit fitting software and put the devices in discovery mode by removing them from the charger or closing the battery door for non-rechargeable models. Click **Connect to ReSound Smart Fit**.
- 2 The hearing aids will appear in the **Select Instruments** section. Use the **Push to Assign** feature to assign to the respective right and left ear placement. Click **Continue** in the lower right of the screen. Select appropriate physical properties and click **Continue** to proceed through connection process.
- 3 Calibrate the hearing aids on the patient's ears to activate DFS Ultra III. Click **Continue** in the lower right of the screen. **Note:** Devices will automatically unmute once calibrated.
- 4 Navigate to **Patient > Profile** screen, customise the Prescribed First Fit selections as appropriate:
 - a Select Patient **Experience Level** based on previous experience with amplification and to maximise acceptance of hearing aids.
 - b Click **Fit Patient** in the lower right, or navigate to **Gain Adjustments**.
- 5 Fine-tune instruments using **Gain Adjustments**:
 - a Add or remove programs from this screen.
 - b Once programming is set, navigate to **Device Controls** to customise the manual push button, demonstrate beeps or pair accessories.
- 6 Click **Save** at bottom right of screen. Devices will reboot and play the start-up melody upon exiting the software.

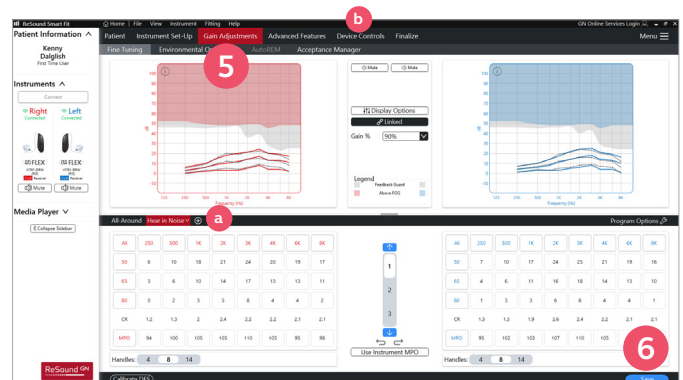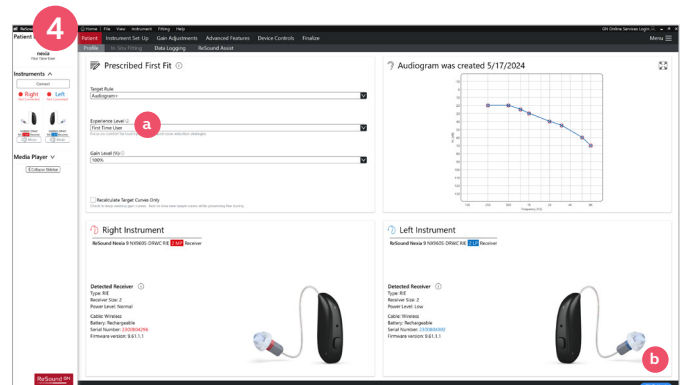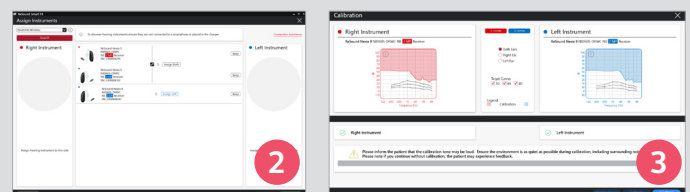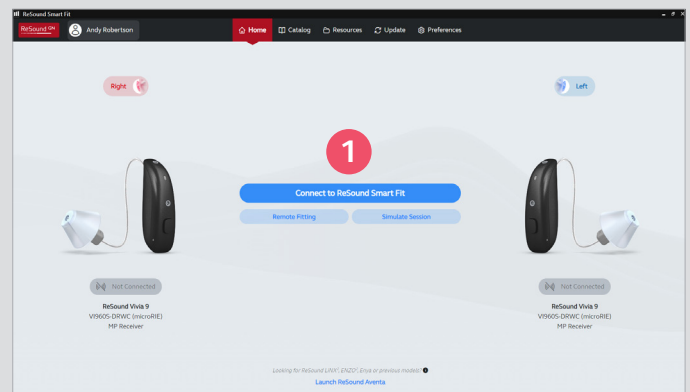

Worldwide Headquarters

GN Hearing A/S
Lautrupbjerg 7
DK-2750 Ballerup
Denmark
Tel.: +45 4575 1111
gnhearing.com

CVR no. 55082715

© 2025 GN Hearing A/S. All rights reserved. ReSound is a trademark of GN Hearing A/S.