



Oticon Opn™ miniRITE is a discreet design with a smart single push button.

Oticon Opn miniRITE T is a discreet style based on the popular miniRITE. It features a telecoil and a convenient double push button for easy operation of the volume and program controls.


OpenSound Navigator™ provides better speech understanding by continuously analyzing the environment, balancing all sound sources and attenuating the dominating noise.

TwinLink™ wireless technology combines binaural communication and 2.4 GHz connectivity in stereo directly to external digital devices with very low power consumption.

Fully programmable with updatable firmware, the Velox™ platform is ready for the future.

Oticon Opn is a Made for iPhone® hearing aid.

Oticon Opn is built on the new Velox platform, providing frequency resolution in 64 channels (Opn 1).

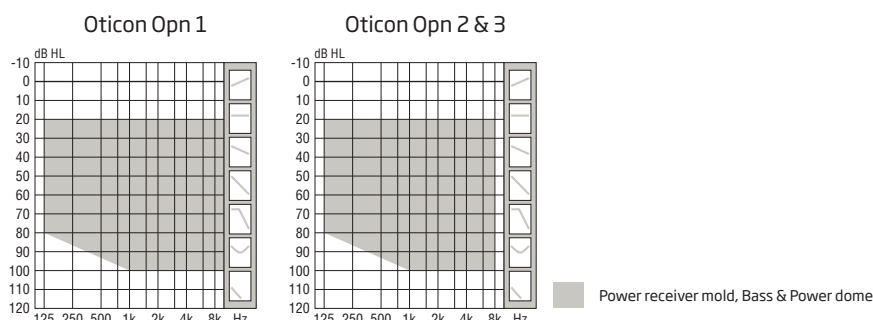

100

		Oticon Opn 1	Oticon Opn 2	Oticon Opn 3
Speech Understanding	OpenSound Navigator™	Level 1	Level 2	Level 3
	- Balancing power effect	100%	50%	50%
	- Max. noise removal	9 dB	5 dB	3 dB
	Speech Guard™ LX	Level 1	Level 2	Level 3
	Spatial Sound™ LX	4 estimators	2 estimators	2 estimators
	Soft Speech Booster LX	•	•	•
	Speech Rescue™ LX	•	•	•
Sound Quality	Clear Dynamics	•	•	-
	Spatial Noise Management	•	•	-
	Fitting Bandwidth*	10 kHz	8 kHz	8 kHz
	Processing Channels	64	48	48
	Bass Boost (streaming)	•	•	•
Listening Comfort	Transient Noise Management	4 configurations	On/Off	On/Off
	Feedback shield LX	•	•	•
	Wind Noise Management	•	•	•
Personalization & Optimizing Fitting	YouMatic™ LX	3 configurations	2 configurations	1 configuration
	Fitting Bands	16	14	12
	Multiple Directionality Options	•	•	•
	Adaptation Management	•	•	•
	Oticon Firmware Updater	•	•	•
	Fitting Formulas	VAC+, NAL-NL1 + 2, DSL v5.0	VAC+, NAL-NL1 + 2, DSL v5.0	VAC+, NAL-NL1 + 2, DSL v5.0
Connecting to the World	Stereo streaming (2.4 GHz)	•	•	•
	Oticon ON app	•	•	•
	ConnectClip	•	•	•
	Remote Control 3.0	•	•	•
	TV Adapter 3.0	•	•	•
	Tinnitus SoundSupport™	•	•	•
Expected battery life, hours**		50-65	50-65	50-65

* Bandwidth accessible for gain adjustments during fitting

** Battery size 312 - IEC PR41.

Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).



IP68

For information on compatibility, please visit www.oticon.com/support/connectivity.


Apple, the Apple logo, iPhone, iPad, and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries.

Technical data Measured according to		Ear Simulator IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010			2CC Coupler ANSI S3.22:2014, IEC 60118-0:2015 and IEC 60318-5:2006		
Oticon Opn miniRITE/miniRITE T		Opn 1	Opn 2	Opn 3	Opn 1	Opn 2	Opn 3
Frequency range Hz		100-8500	100-7500	100-7500	100-8000	100-7500	100-7500
OSPL90	Peak	132 dB SPL			122 dB SPL		
	1600 Hz	130 dB SPL			121 dB SPL		
	HFA-OSPL90	127 dB SPL			118 dB SPL		
Full-on gain*	Peak	66 dB			57 dB		
	1600 Hz	56 dB			48 dB		
	HFA-FOG	59 dB			51 dB		
Reference test gain		49 dB			42 dB		
Telecoil output (1600 Hz) (miniRITE T)	1 mA/m field	86 dB SPL			-		
	10 mA/m field	106 dB SPL			-		
	SPLITS L/R	-			103/103 dB SPL		
Total harmonic distortion (Input 70 dB SPL)	500 Hz	< 7 %			< 2 %		
	800 Hz	< 4 %			< 2 %		
	1600 Hz	< 2 %			< 2 %		
Equivalent input noise level	Omni	23 dB SPL			19 dB SPL		
	Dir	32 dB SPL			30 dB SPL		
Battery consumption**	Typical	1.5 mA			1.7 mA		
	Quiescent	1.5 mA			1.5 mA		
Battery life, artificial measurement, hours***		115			105		

* Measured with the gain control of the hearing aid set to its full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0+A1:1994 but without influence of feedback.

** Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.

*** Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.

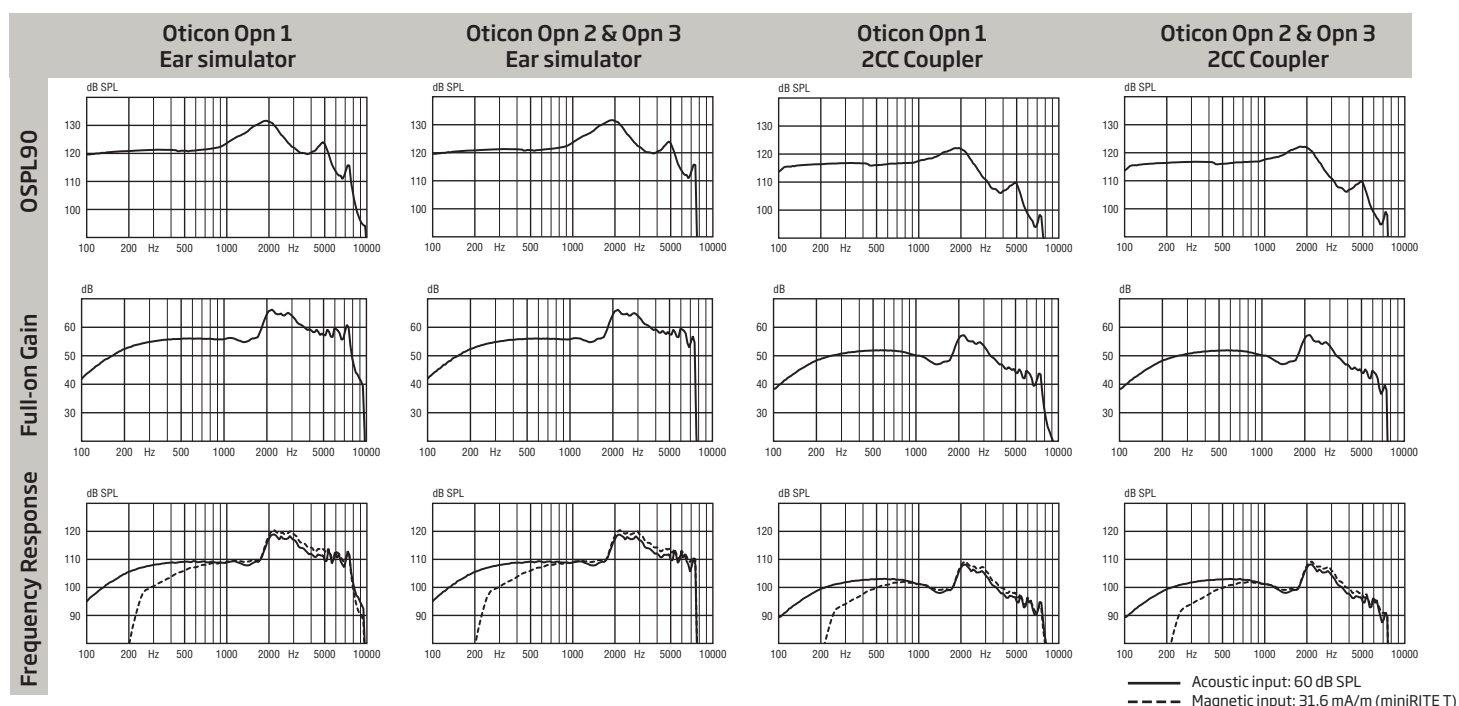

Technical information: Omnidirectional mode is used unless otherwise stated.

Operating conditions

Temperature: +1°C to +40°C
(34°F to 104°F)

Humidity: 5% - 93% relative
humidity, non-condensing
Atmospheric pressure:
700 hPa to 1060 hPa

Storage and transportation conditions

Temperature and humidity should not exceed the
below limits for extended periods during
transportation and storage.

Temperature: -25°C to +60°C (-13°F to 140°F)
Humidity: 5% - 93% relative humidity, non-condensing
Atmospheric pressure: 700 hPa to 1060 hPa

Warning to the hearing aid dispenser

The maximum output capability of the hearing
aid may exceed 132 dB SPL (IEC 711). Special care
should be exercised in selecting and fitting the
instrument as there may be risk of impairing the
remaining hearing of the hearing instrument user.



SBO Hearing A/S
Kongebakken 9
DK-2765 Smørum
Denmark

Headquarters
Oticon A/S
Kongebakken 9
DK-2765 Smørum, Denmark

www.oticon.com

Oticon is part of the Demant Group.

oticon
life-changing technology