



	Oticon Opn Play 1	Oticon Opn Play 2	
Speech Understanding	OpenSound Navigator™	Level 1	Level 3
	- Balancing power effect	100%	50%
	- Max. noise removal	9 dB	3 dB
	OpenSound Optimizer™	•	•
	Speech Guard™ LX	Level 1	Level 3
	Spatial Sound™ LX	4 estimators	2 estimators
	Speech Rescue™ LX	•	•
Sound Quality	Clear Dynamics	•	-
	Spatial Noise Management	•	-
	Fitting Bandwidth*	10 KHz	8 KHz
	Processing Channels	64	48
	Bass Boost (streaming)	•	•
Listening Comfort	Transient Noise Management	4 configurations	On/Off
	Feedback shield LX	•	•
	Wind Noise Management	•	•
Optimizing Fitting	YouMatic™ LX	3 configurations	1 configuration
	Fitting Bands	16	12
	REM AutoFit	Verifit®LINK, IMC 2**	Verifit®LINK, IMC 2**
	Pediatric Fitting Mode	•	•
	DSL Fitting Range ***	•	•
	Fitting Formulas	DSL v5.0, NAL-NL1 + 2, VAC+	DSL v5.0, NAL-NL1 + 2, VAC+
Designed for children	Hypo Allergenic	•	•
	IP Rating	IP68	IP68
	Nano Coating	•	•
	Color Options	12	12
	Integrated 2.4 GHz Receiver	•	•

* Bandwidth accessible for gain adjustments during fitting

** Inter Module Communication 2

*** Available in this Technical Data sheet and Opn Play Product Guide

Operating conditions

Temperature: +1°C to +40°C (34°F to 104°F)

Humidity: 5% to 93% relative humidity, non-condensing

Atmospheric pressure: 700 hPa to 1060 hPa

Storage and transportation conditions

Temperature and humidity shall not exceed the below limits for extended periods during transportation and storage.

Temperature: -25°C to +60°C (-13°F to 140°F)

Humidity: 5% to 93% relative humidity, non-condensing

Atmospheric pressure: 700 hPa to 1060 hPa

Oticon Opn Play™ miniRITE is a stylish and discreet hearing aid that sits snugly behind the ear. Connectivity and remote microphone access are available through integrated 2.4 GHz technology.

OpenSound Navigator™ gives children 360° access to the full soundscape across simple and complex listening environments, constantly maximizing learning opportunities.

OpenSound Optimizer™ proactively identifies and prevents feedback before it occurs, preserving audibility and increasing fitting flexibility.

TwinLink™ wireless technology combines binaural communication and 2.4 GHz connectivity with stereo streaming directly from digital devices.

Oticon Opn Play is powered by Oticon's Velox S™ platform which offers market-leading speed and resolution and takes pediatric hearing care to a new level.

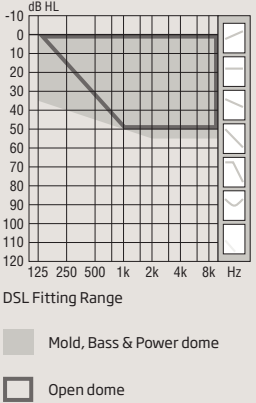
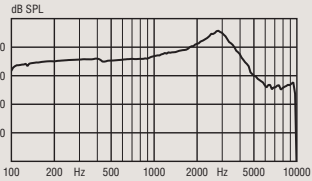
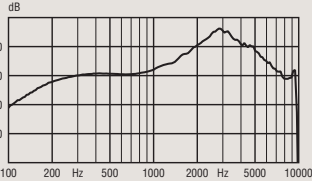
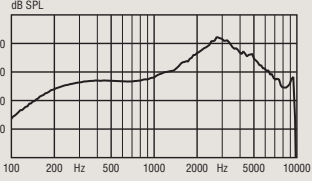
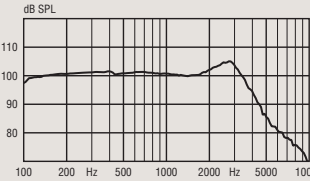
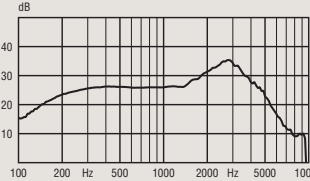
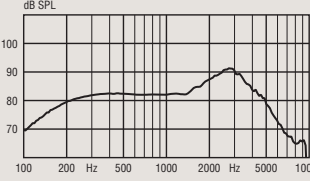
Apple, the Apple logo, iPhone, iPad, and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries.



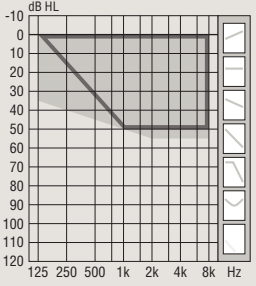
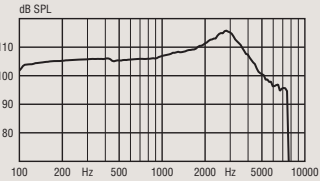
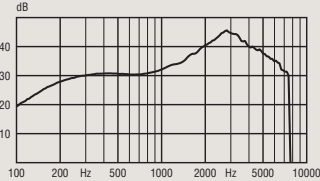
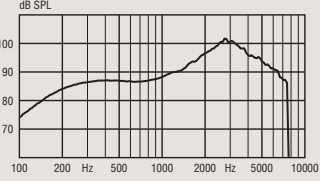
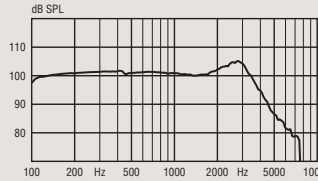
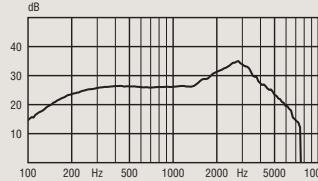
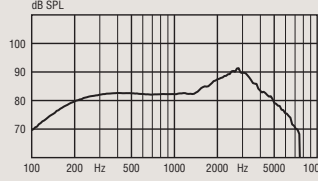
IP68



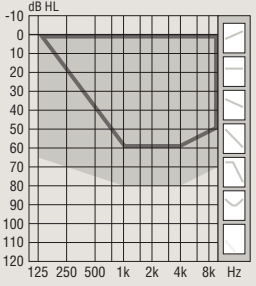
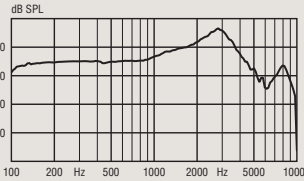
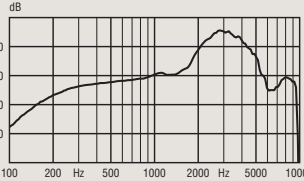
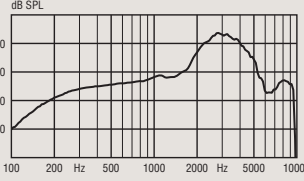
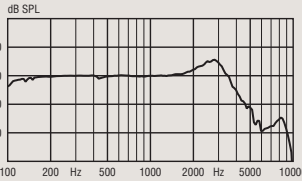
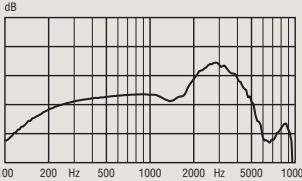
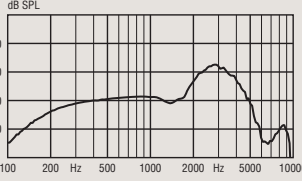
For information on compatibility, please visit www.oticon.com/support/compatibility

Technical data		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 60  Technical information Omnidirectional mode is used unless otherwise stated.		OSPL90  Full-on Gain  Frequency Response 	OSPL90  Full-on Gain  Frequency Response 
OSPL90	Peak 1600 Hz HFA-OSPL90	116 dB SPL 109 dB SPL 110 dB SPL	105 dB SPL 100 dB SPL 102 dB SPL
Full-on gain ¹	Peak 1600 Hz HFA-FOG	46 dB 37 dB 38 dB	35 dB 29 dB 30 dB
Reference test gain		30 dB	26 dB
Frequency range		110-9700 Hz	100-9200 Hz
Telecoil output (1600 Hz)	1 mA/m field 10 mA/m field SPLITS L/R	- - -	- - -
Total harmonic distortion (Input 70 dB SPL)	500 Hz 800 Hz 1600 Hz	< 2 % < 3 % < 2 %	< 2 % < 2 % < 2 %
Equivalent input noise level	Omni Dir	21 dB SPL 28 dB SPL	18 dB SPL 27 dB SPL
Battery consumption ²	Typical Quiescent	1.5 mA 1.5 mA	1.6 mA 1.5 mA
Battery life, artificial measurement, hours ³		120	115
Expected battery life, hours (battery size 312 - IEC PR41) ⁴		60-65	

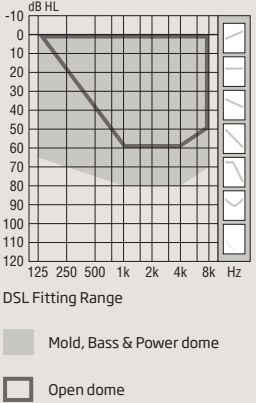
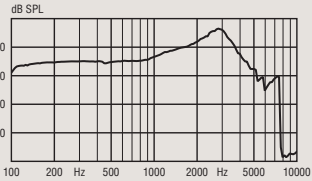
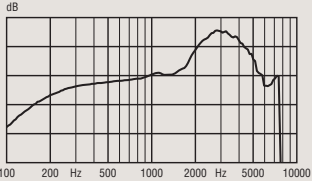
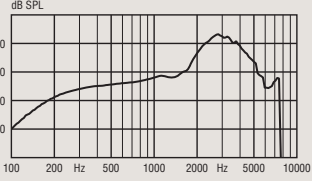
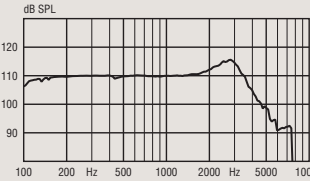
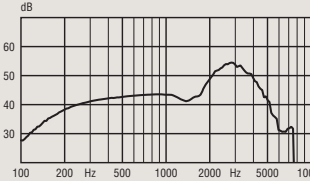
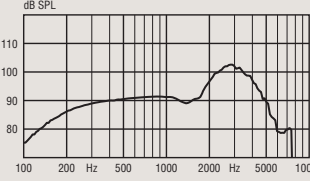
- 1) Measured with the gain control of the hearing aid set to its full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0+A1:1994 but without influence of feedback.
- 2) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.
- 3) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
- 4) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

Technical data		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 60  DSL Fitting Range ■ Mold, Bass & Power dome □ Open dome Technical information Omnidirectional mode is used unless otherwise stated.		OSPL90  Full-on Gain  Frequency Response 	OSPL90  Full-on Gain  Frequency Response 
OSPL90	Peak 1600 Hz HFA-OSPL90	116 dB SPL 109 dB SPL 110 dB SPL	105 dB SPL 100 dB SPL 102 dB SPL
Full-on gain ¹	Peak 1600 Hz HFA-FOG	46 dB 37 dB 38 dB	35 dB 29 dB 30 dB
Reference test gain		30 dB	26 dB
Frequency range		110-7500 Hz	100-7500 Hz
Telecoil output (1600 Hz)	1 mA/m field 10 mA/m field SPLITS L/R	- - -	- - -
Total harmonic distortion (Input 70 dB SPL)	500 Hz 800 Hz 1600 Hz	< 2 % < 3 % < 2 %	< 2 % < 2 % < 2 %
Equivalent input noise level	Omni Dir	22 dB SPL 30 dB SPL	19 dB SPL 28 dB SPL
Battery consumption ²	Typical Quiescent	1.5 mA 1.5 mA	1.6 mA 1.5 mA
Battery life, artificial measurement, hours ³		120	115
Expected battery life, hours (battery size 312 - IEC PR41) ⁴		60-65	

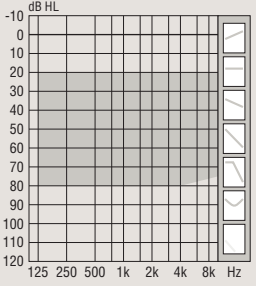
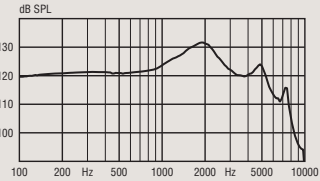
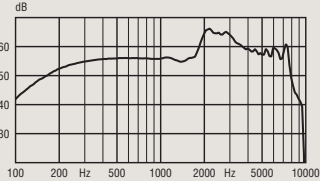
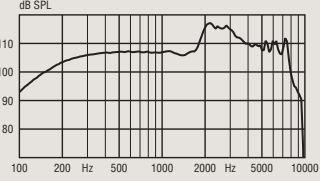
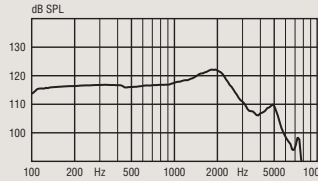
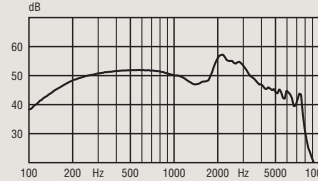
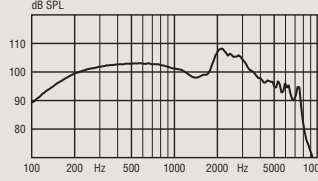
- 1) Measured with the gain control of the hearing aid set to its full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0+A1:1994 but without influence of feedback.
- 2) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.
- 3) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
- 4) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

Technical data		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 DSL Fitting Range ■ Mold, Bass & Power dome □ Open dome Technical information Omnidirectional mode is used unless otherwise stated.		OSPL90  Full-on Gain  Frequency Response 	OSPL90  Full-on Gain  Frequency Response 
OSPL90	Peak 1600 Hz HFA-OSPL90	127 dB SPL 120 dB SPL 121 dB SPL	116 dB SPL 111 dB SPL 112 dB SPL
Full-on gain ¹	Peak 1600 Hz HFA-FOG	66 dB 52 dB 55 dB	54 dB 43 dB 47 dB
Reference test gain		45 dB	34 dB
Frequency range		120-9500 Hz	100-8500 Hz
Telecoil output (1600 Hz)	1 mA/m field 10 mA/m field SPLITS L/R	- - -	- - -
Total harmonic distortion (Input 70 dB SPL)	500 Hz 800 Hz 1600 Hz	< 2 % < 3 % < 2 %	< 2 % < 2 % < 2 %
Equivalent input noise level	Omni Dir	25 dB SPL 32 dB SPL	20 dB SPL 29 dB SPL
Battery consumption ²	Typical Quiescent	1.6 mA 1.5 mA	1.7 mA 1.5 mA
Battery life, artificial measurement, hours ³		110	105
Expected battery life, hours (battery size 312 - IEC PR41) ⁴		55-65	

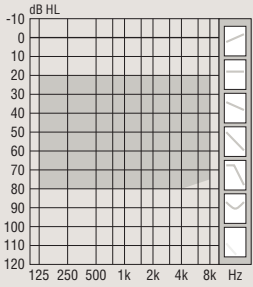
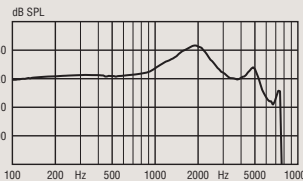
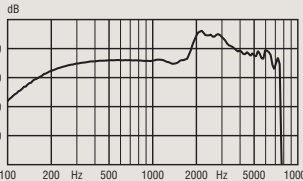
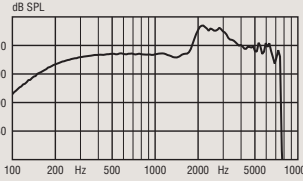
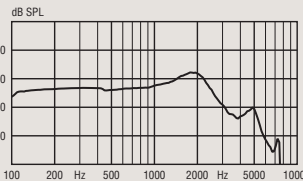
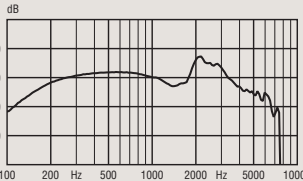
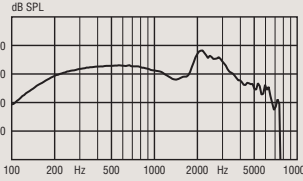
- 1) Measured with the gain control of the hearing aid set to its full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0+A1:1994 but without influence of feedback.
- 2) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.
- 3) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
- 4) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

Technical data		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 85  Technical information Omnidirectional mode is used unless otherwise stated.		OSPL90  Full-on Gain  Frequency Response 	OSPL90  Full-on Gain  Frequency Response 
OSPL90		Peak 1600 Hz HFA-OSPL90	127 dB SPL 120 dB SPL 121 dB SPL
Full-on gain ¹		Peak 1600 Hz HFA-FOG	66 dB 52 dB 55 dB
Reference test gain			45 dB
Frequency range			120-7500 Hz
Telecoil output (1600 Hz)		1 mA/m field 10 mA/m field SPLITS L/R	- - -
Total harmonic distortion (Input 70 dB SPL)		500 Hz 800 Hz 1600 Hz	< 2 % < 3 % < 2 %
Equivalent input noise level		Omni Dir	26 dB SPL 33 dB SPL
Battery consumption ²		Typical Quiescent	1.6 mA 1.5 mA
Battery life, artificial measurement, hours ³			110
Expected battery life, hours (battery size 312 - IEC PR41) ⁴			55-65

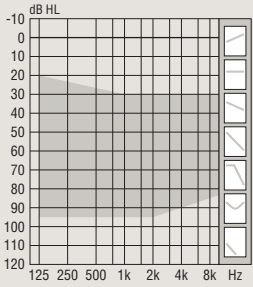
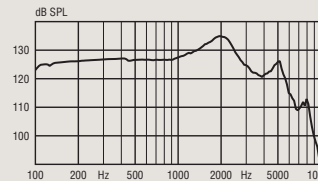
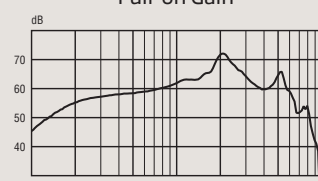
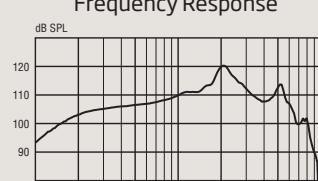
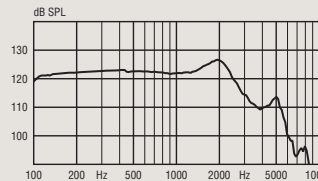
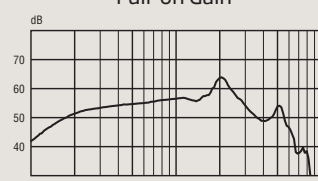
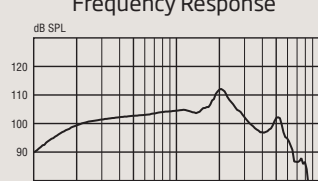
- 1) Measured with the gain control of the hearing aid set to its full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0+A1:1994 but without influence of feedback.
- 2) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.
- 3) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
- 4) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

Technical data		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 DSL Fitting Range  Power receiver mold, Bass & Power dome		OSPL90  Full-on Gain  Frequency Response 	OSPL90  Full-on Gain  Frequency Response 
OSPL90	Peak 1600 Hz HFA-OSPL90	132 dB SPL 130 dB SPL 127 dB SPL	122 dB SPL 121 dB SPL 118 dB SPL
Full-on gain ¹	Peak 1600 Hz HFA-FOG	66 dB 56 dB 59 dB	57 dB 48 dB 51 dB
Reference test gain		49 dB	42 dB
Frequency range		100-8500 Hz	100-8000 Hz
Telecoil output (1600 Hz)	1 mA/m field 10 mA/m field SPLITS L/R	- - -	- - -
Total harmonic distortion (Input 70 dB SPL)	500 Hz 800 Hz 1600 Hz	< 7 % < 4 % < 2 %	< 2 % < 2 % < 2 %
Equivalent input noise level	Omni Dir	23 dB SPL 32 dB SPL	19 dB SPL 30 dB SPL
Battery consumption ²	Typical Quiescent	1.5 mA 1.5 mA	1.7 mA 1.5 mA
Battery life, artificial measurement, hours ³		115	105
Expected battery life, hours (battery size 312 - IEC PR41) ⁴		50-65	

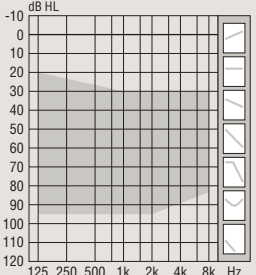

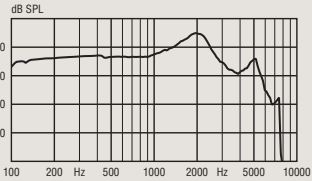
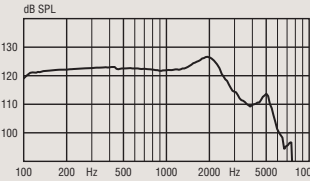
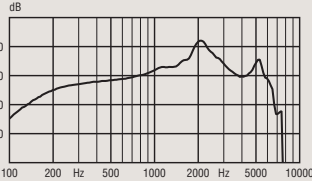
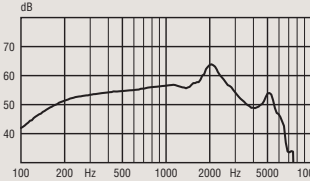
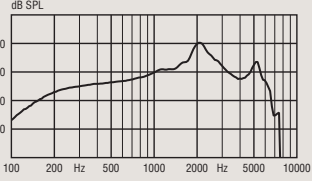
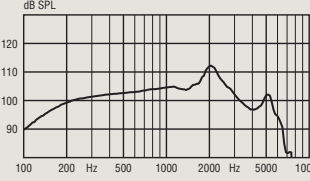
- 1) Measured with the gain control of the hearing aid set to its full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0+A1:1994 but without influence of feedback.
- 2) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.
- 3) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
- 4) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

Technical data		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 DSL Fitting Range Power receiver mold, Bass & Power dome  Technical information Omnidirectional mode is used unless otherwise stated. Warning to the hearing care provider The maximum output capability of the hearing aid may exceed 132 dB SPL (IEC 711). Special care should be exercised in selecting and fitting the instrument as there may be risk of impairing the remaining hearing of the hearing aid user.		OSPL90  Full-on Gain  Frequency Response 	OSPL90  Full-on Gain  Frequency Response 
OSPL90	Peak 1600 Hz HFA-OSPL90	132 dB SPL 130 dB SPL 127 dB SPL	122 dB SPL 121 dB SPL 118 dB SPL
Full-on gain ¹	Peak 1600 Hz HFA-FOG	66 dB 56 dB 59 dB	57 dB 48 dB 51 dB
Reference test gain		49 dB	42 dB
Frequency range		100-7500 Hz	100-7500 Hz
Telecoil output (1600 Hz)	1 mA/m field 10 mA/m field SPLITS L/R	- - -	- - -
Total harmonic distortion (Input 70 dB SPL)	500 Hz 800 Hz 1600 Hz	< 7 % < 4 % < 2 %	< 2 % < 2 % < 2 %
Equivalent input noise level	Omni Dir	23 dB SPL 32 dB SPL	19 dB SPL 30 dB SPL
Battery consumption ²	Typical Quiescent	1.5 mA 1.5 mA	1.7 mA 1.5 mA
Battery life, artificial measurement, hours ³		115	105
Expected battery life, hours (battery size 312 - IEC PR41) ⁴		50-65	

- 1) Measured with the gain control of the hearing aid set to its full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0+A1:1994 but without influence of feedback.
- 2) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.
- 3) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
- 4) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

Technical data		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 105 DSL Fitting Range Power receiver mold		OSPL90  Full-on Gain  Frequency Response 	OSPL90  Full-on Gain  Frequency Response 
OSPL90	Peak 1600 Hz HFA-OSPL90	135 dB SPL 132 dB SPL 130 dB SPL	127 dB SPL 125 dB SPL 122 dB SPL
Full-on gain ¹	Peak 1600 Hz HFA-FOG	72 dB 65 dB 65 dB	64 dB 57 dB 57 dB
Reference test gain		58 dB	46 dB
Frequency range		100-8200 Hz	100-7800 Hz
Telecoil output (1600 Hz)	1 mA/m field 10 mA/m field SPLITS L/R	- - -	- - -
Total harmonic distortion (Input 70 dB SPL)	500 Hz 800 Hz 1600 Hz	<2 % <2 % <3 %	<2 % <2 % <2 %
Equivalent input noise level	Omni Dir	18 dB SPL 28 dB SPL	18 dB SPL 29 dB SPL
Battery consumption ²	Typical Quiescent	1.6 mA 1.5 mA	1.7 mA 1.5 mA
Battery life, artificial measurement, hours ³		110	105
Expected battery life, hours (battery size 312 - IEC PR41) ⁴		45-65	

- 1) Measured with the gain control of the hearing aid set to its full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0+A1:1994 but without influence of feedback.
- 2) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.
- 3) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
- 4) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

Technical data		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 DSL Fitting Range Power receiver mold		OSPL90 	OSPL90 
		Full-on Gain 	Full-on Gain 
		Frequency Response 	Frequency Response 
OSPL90	Peak	135 dB SPL	127 dB SPL
	1600 Hz	132 dB SPL	125 dB SPL
	HFA-OSPL90	130 dB SPL	122 dB SPL
Full-on gain ¹	Peak	72 dB	64 dB
	1600 Hz	65 dB	57 dB
	HFA-FOG	65 dB	57 dB
Reference test gain		58 dB	46 dB
Frequency range		100-7500 Hz	100-6500 Hz
Telecoil output (1600 Hz)	1 mA/m field	-	-
	10 mA/m field	-	-
	SPLITS L/R	-	-
Total harmonic distortion (Input 70 dB SPL)	500 Hz	< 2 %	< 2 %
	800 Hz	< 2 %	< 2 %
	1600 Hz	< 3 %	< 2 %
Equivalent input noise level	Omni	18 dB SPL	18 dB SPL
	Dir	28 dB SPL	29 dB SPL
Battery consumption ²	Typical	1.6 mA	1.7 mA
	Quiescent	1.5 mA	1.5 mA
Battery life, artificial measurement, hours ³		110	105
Expected battery life, hours (battery size 312 - IEC PR41) ⁴		45-65	

- Measured with the gain control of the hearing aid set to its full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0+A1:1994 but without influence of feedback.
- Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.
- Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
- Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

Headquarters
Oticon A/S
Kongebakken 9
DK-2765 Smørum
Denmark



SBO Hearing A/S
Kongebakken 9
DK-2765 Smørum
Denmark

237500US / 2022.03.30 / v1