

Technical data sheet

BTE PP

105



	Oticon Opn Play 1	Oticon Opn Play 2	
Speech Understanding	OpenSound Navigator™	Level 1	Level 3
	- Balancing power effect	100%	50%
	- Max. noise removal	9 dB	3 dB
	OpenSound Optimizer™	•	•
	Speech Guard™ LX	Level 1	Level 3
	Spatial Sound™ LX	4 estimators	2 estimators
	Speech Rescue™ LX	•	•
Sound Quality	Clear Dynamics	•	-
	Spatial Noise Management	•	-
	Fitting Bandwidth*	10 KHz	8 KHz
	Processing Channels	64	48
	Bass Boost (streaming)	•	•
Listening Comfort	Transient Noise Management	4 configurations	On/Off
	Feedback shield LX	•	•
	Wind Noise Management	•	•
Optimizing Fitting	YouMatic™ LX	3 configurations	1 configuration
	Fitting Bands	16	12
	REM AutoFit	Verifit®LINK, IMC 2**	Verifit®LINK, IMC 2**
	Pediatric Fitting Mode	•	•
	DSL Fitting Range ***	•	•
	Fitting Formulas	DSL v5.0, NAL-NL1 + 2, VAC+	DSL v5.0, NAL-NL1 + 2, VAC+
Designed for children	LED	•	•
	Tamper Resistant Battery Drawer	•	•
	Hypo Allergenic	•	•
	IP Rating	IP68	IP68
	Nano Coating	•	•
	Color Options	12	12
	Integrated 2.4 GHz Receiver	•	•
	DAI/FM Compatibility	•	•

* Bandwidth accessible for gain adjustments during fitting

** Inter Module Communication 2

*** Available in this Technical Data sheet and Opn Play Product Guide

Operating conditions

Temperature: +34°F to +104°F

Relative humidity: 5% to 93%, non-condensing

Storage and transportation conditions

Temperature and humidity should not exceed the following limits for extended periods during transportation and storage.

Temperature: -13°F to +140°F

Relative humidity: 5% to 93%, non-condensing

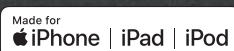
Oticon Opn Play™ BTE PP has a compact design and offers both hook and thin tube. Features telecoil, double push button with optional LED status and FM compatibility.

OpenSound Navigator™ helps users to select and understand speech in all types of environments by balancing the sound sources and attenuating noise.

OpenSound Optimizer™ improves users listening experience and comfort by blocking feedback and securing the targeted amplification of sound sources.

TwinLink™ wireless technology combines binaural communication and 2.4 GHz connectivity with stereo streaming directly from digital devices.

Oticon Opn Play is built on the powerful Velox S™ platform which has a programmable firmware architecture, supporting future performance updates.



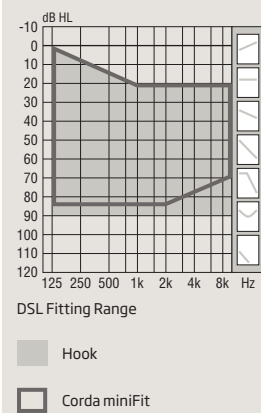
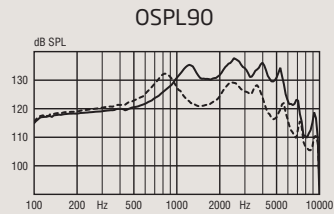
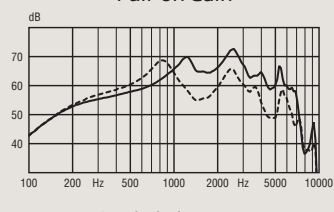
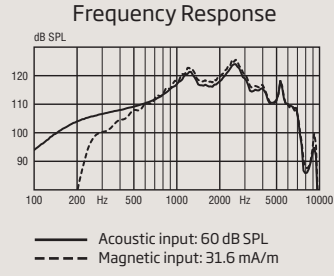
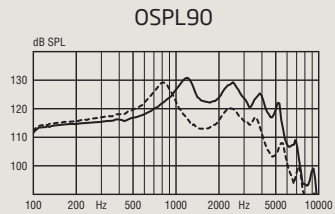
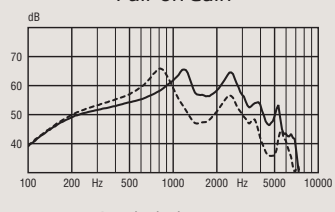
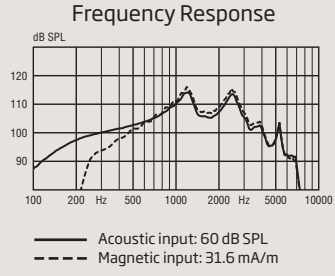
IP68

oticon
PEDIATRICS

For information on compatibility, please visit www.oticon.com/support

Oticon Opn Play 1

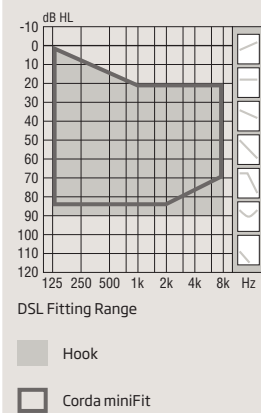
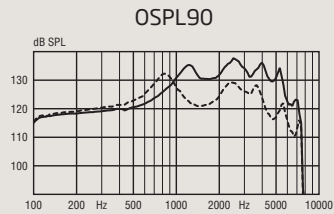
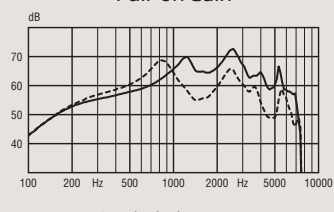
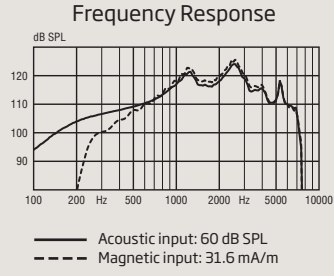
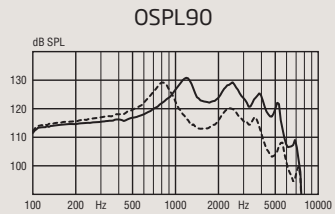
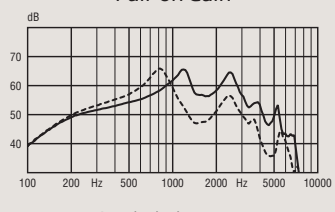
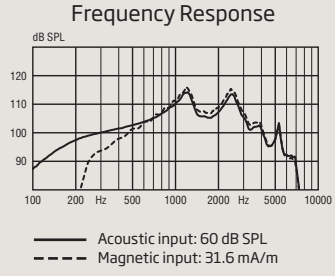
BTE PP

Technical data		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 105 		  	  
OSPL90	Peak 1600 Hz HFA-OSPL90	138 (132 ¹) dB SPL 130 (121 ¹) dB SPL 133 (126 ¹) dB SPL	131 (129 ¹) dB SPL 123 (113 ¹) dB SPL 126 (118 ¹) dB SPL
Full-on gain ²	Peak 1600 Hz HFA-FOG	73 (69 ¹) dB 65 (56 ¹) dB 68 (62 ¹) dB	66 (66 ¹) dB 57 (47 ¹) dB 61 (54 ¹) dB
Reference test gain		57 dB	50 dB
Frequency range		150-7300 Hz	120-7000 Hz
Telecoil output (1600 Hz)	1 mA/m field 10 mA/m field SPLITS L/R	97 dB SPL 117 dB SPL -	- - 109/109 dB SPL
Total harmonic distortion (Input 70 dB SPL)	500 Hz 800 Hz 1600 Hz	7 % 5 % < 2 %	3 % < 2 % < 2 %
Equivalent input noise level	Omni Dir	17 dB SPL 29 dB SPL	14 dB SPL 27 dB SPL
Battery consumption ³	Typical Quiescent	1.8 mA 1.6 mA	1.9 mA 1.6 m
Battery life, artificial measurement, hours ⁴		175	160
Expected battery life, hours (battery size 13 - IEC PR48) ⁵ IRIL (IEC 60118-13:2011)		80-105 700/1400/2000 MHz: 18/20/40 dB SPL	

1) For instruments fitted with Corda miniFit Power.
2) Measured with the gain control of the hearing aid set to its full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0+A1:1994 but without influence of feedback.
3) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.
4) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
5) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

Oticon Opn Play 2

BTE PP

Technical data		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 105 		  	  
OSPL90	Peak 1600 Hz HFA-OSPL90	138 (132 ¹) dB SPL 130 (121 ¹) dB SPL 133 (126 ¹) dB SPL	131 (129 ¹) dB SPL 123 (113 ¹) dB SPL 126 (118 ¹) dB SPL
Full-on gain ²	Peak 1600 Hz HFA-FOG	73 (69 ¹) dB 65 (56 ¹) dB 68 (62 ¹) dB	66 (66 ¹) dB 57 (47 ¹) dB 61 (54 ¹) dB
Reference test gain		57 dB	50 dB
Frequency range		150-7300 Hz	120-7000 Hz
Telecoil output (1600 Hz)	1 mA/m field 10 mA/m field SPLITS L/R	97 dB SPL 117 dB SPL -	- - 109/109 dB SPL
Total harmonic distortion (Input 70 dB SPL)	500 Hz 800 Hz 1600 Hz	7 % 5 % < 2 %	3 % < 2 % < 2 %
Equivalent input noise level	Omni Dir	17 dB SPL 29 dB SPL	14 dB SPL 27 dB SPL
Battery consumption ³	Typical Quiescent	1.8 mA 1.6 mA	1.9 mA 1.6 m
Battery life, artificial measurement, hours ⁴		175	160
Expected battery life, hours (battery size 13 - IEC PR48) ⁵ IRIL (IEC 60118-13:2011)		80-105 700/1400/2000 MHz: 18/20/40 dB SPL	

1) For instruments fitted with Corda miniFit Power.
2) Measured with the gain control of the hearing aid set to its full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0+A1:1994 but without influence of feedback.
3) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.
4) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
5) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).



Oticon A/S
Kongebakken 9
2765 Smørum
Denmark
+45 3917 7100

198298US/2018.11.24



www.oticon.com

oticon
PEDIATRICS