



	Play PX 1	Play PX 2	
Speech Understanding	MoreSound Intelligence™	Level 1	Level 3
	- Environment configuration	5 Options	3 Options
	- Virtual Outer Ear	3 Configurations	1 Configuration
	- Spatial Balancer	100%	60%
	- Neural Noise Suppression, Difficult / Easy	10 dB / 4 dB	6 dB / 0 dB
	- Sound Enhancer	3 Configurations	1 Configuration
	MoreSound Amplifier™	•	•
	Feedback Prevention	MoreSound Optimizer™ & Feedback shield	MoreSound Optimizer™ & Feedback shield
	Spatial Sound™	4 Estimators	2 Estimators
	Soft Speech Booster	•	•
Frequency lowering	Speech Rescue™	Speech Rescue™	
Sound Quality	Clear Dynamics	•	-
	Better-Ear Priority	•	-
	Fitting Bandwidth*	10 kHz	8 kHz
	Bass Boost (streaming)	•	•
	Processing Channels	64	48
Listening Comfort	Transient Noise Management	4 configurations	3 configurations
	Wind Noise Management	•	•
Optimizing Fitting	Fitting Bands	24	18
	REM Autofit	Verifit®LINK, IMC 2**	Verifit®LINK, IMC 2**
	Pediatric Fitting Mode	•	•
	DSL Fitting Range***	•	•
	Fitting Formulas	DSL v5.0, NAL-NL 1/ NAL-NL 2, VAC+	DSL v5.0, NAL-NL 1/ NAL-NL 2, VAC+
Designed for children	Tamper Resistant Battery Drawer	•	•
	LED	•	•
	Biologically safe	•	•
	Nano coating	•	•
	Color options	12	12
	Hands-free communication****	•	•
	Direct streaming*****	•	•
	Edumic	•	•
	Oticon ON app	•	•

* Bandwidth accessible for gain adjustments during fitting
 ** Inter Module Communication 2

*** Available in this Technical Data sheet and Oticon Play PX Product Guide

**** Available for Oticon Play PX from FW 1.1 with selected iPhone® models

***** From iPhone®, iPad®, iPod touch®, and selected Android™ devices

Operating Conditions

Temperature: +1°C to +40°C (34°F to 104°F)

Humidity: 5% to 93% relative humidity, non-condensing

Atmospheric pressure: 700 hPa to 1060 hPa

Storage and transportation conditions

Temperature and humidity should not exceed the below limits for extended periods during transportation and storage.

Transportation

Temperature: -25°C to +60°C (-13°F to 140°F)

Humidity: 5% to 93% relative humidity, non-condensing

Atmospheric pressure: 700 hPa to 1060 hPa

Storage

Temperature: -25°C to +60°C (-13°F to 140°F)

Humidity: 5% to 93% relative humidity, non-condensing

Atmospheric pressure: 700 hPa to 1060 hPa

Oticon Play PX miniBTE T is small in size and fits most ears. It comes with an LED-light for easy handling. The style features telecoil and a single push-button, and it is powered by a disposable zinc-air battery. It is a Made for iPhone® hearing aid and compatible with the new Android™ protocol for Audio Streaming for Hearing Aids (ASHA) - making it possible to stream directly from iPhone®, iPad®, iPod touch® and select Android devices.

MoreSound Intelligence™ creates a more precise and natural representation of individual sounds with clearer and more distinct contrasts providing access to all relevant sounds.

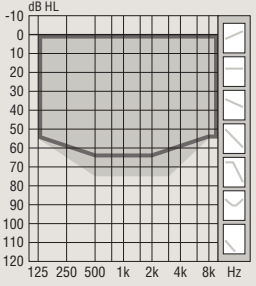
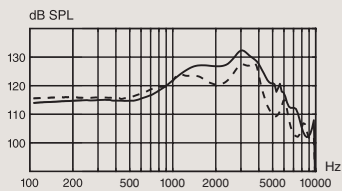
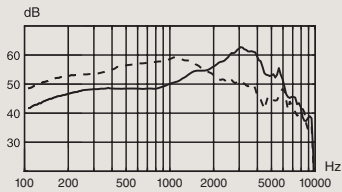
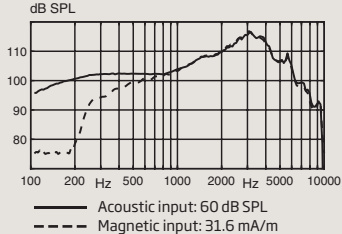
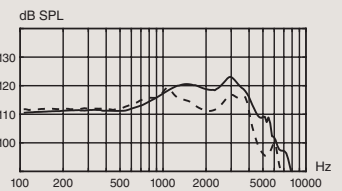
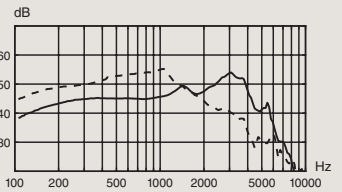
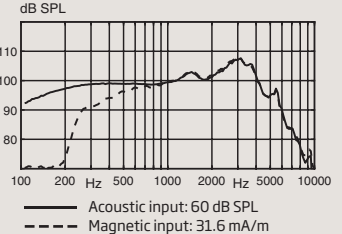
MoreSound Amplifier™ analyzes details in sound, and optimally amplifies them for the brain to have access to relevant information.

Oticon Play PX is built on the innovative Polaris™ platform, which uses a Deep Neural Network to rapidly and optimally manage incoming sounds based on individual needs.

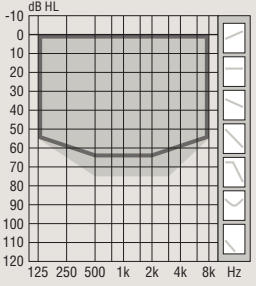
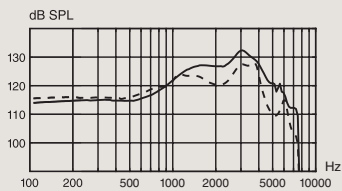
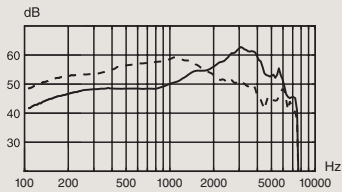
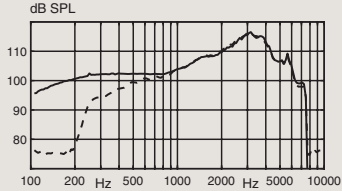
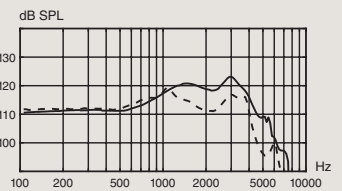
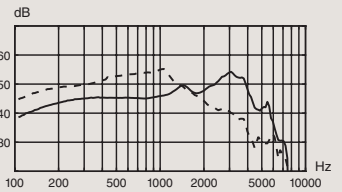
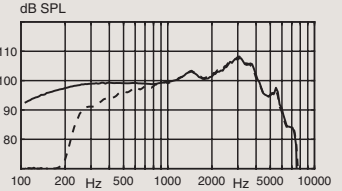
Apple, the Apple logo, iPhone®, iPad®, and iPod touch® are trademarks of Apple Inc., registered in the U.S. and other countries.



For information on compatibility, please visit www.oticon.com/support/compatibility

		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 DSL Fitting Range Hook Corda miniFit		OSPL90  Full-on gain  Frequency response 	OSPL90  Full-on gain  Frequency response 
OSPL90		Peak 132 (128 ¹) dB SPL 1600 Hz 127 (123 ¹) dB SPL HFA-OSPL90 126 (122 ¹) dB SPL	Peak 123 (119 ¹) dB SPL 1600 Hz 120 (114 ¹) dB SPL HFA-OSPL90 119 (115 ¹) dB SPL
Full-on gain ²		Peak 63 (59 ¹) dB 1600 Hz 55 (56 ¹) dB HFA-FOG 55 (55 ¹) dB	Peak 54 (55 ¹) dB 1600 Hz 48 (48 ¹) dB HFA-FOG 48 (48 ¹) dB
Reference test gain		48 dB	42 dB
Frequency range		100-9500 Hz	100-7300 Hz
Telecoil output (1600 Hz)		1 mA/m field 86 dB SPL 10 mA/m field 106 dB SPL SPLITS L/R -	-
Total harmonic distortion (Input 70 dB SPL)		500 Hz < 4 % 800 Hz < 4 % 1600 Hz < 2 %	100/100 dB SPL < 4 % < 3 % < 2 %
Equivalent input noise level		Omni 18 dB SPL Dir 28 dB SPL	17 dB SPL 29 dB SPL
Battery consumption ³		Typical 1.9 mA Quiescent 1.9 mA	2.0 mA 1.9 mA
Battery life, artificial measurement, hours ⁴		95	90
Expected battery life, hours (battery size 312 - IEC PR41) ⁵		50-55	

- 1) For instruments fitted with Corda miniFit Power
- 2) Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0:1983+A1:1994 but without influence of feedback.
- 3) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of a minimum of 3 minutes.
- 4) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
- 5) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 DSL Fitting Range Hook Corda miniFit		OSPL90  Full-on gain  Frequency response 	OSPL90  Full-on gain  Frequency response 
OSPL90		Peak 132 (128 ¹) dB SPL 1600 Hz 127 (123 ¹) dB SPL HFA-OSPL90 126 (122 ¹) dB SPL	Peak 123 (119 ¹) dB SPL 1600 Hz 121 (114 ¹) dB SPL HFA-OSPL90 119 (115 ¹) dB SPL
Full-on gain ²		Peak 63 (59 ¹) dB 1600 Hz 55 (56 ¹) dB HFA-FOG 55 (55 ¹) dB	Peak 54 (55 ¹) dB 1600 Hz 48 (48 ¹) dB HFA-FOG 48 (48 ¹) dB
Reference test gain		48 dB	42 dB
Frequency range		100-7500 Hz	100-7300 Hz
Telecoil output (1600 Hz)		1 mA/m field 86 dB SPL 10 mA/m field 106 dB SPL	-
SPLITS L/R		-	100/100 dB SPL
Total harmonic distortion (Input 70 dB SPL)		500 Hz < 4 % 800 Hz < 4 % 1600 Hz < 2 %	< 4 % < 3 % < 2 %
Equivalent input noise level		Omni 18 dB SPL Dir 28 dB SPL	17 dB SPL 29 dB SPL
Battery consumption ³		Typical 1.9 mA Quiescent 1.9 mA	2.0 mA 1.9 mA
Battery life, artificial measurement, hours ⁴		95	90
Expected battery life, hours (battery size 312 - IEC PR41) ⁵		50-55	

- 1) For instruments fitted with Corda miniFit Power
- 2) Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0:1983+A1:1994 but without influence of feedback.
- 3) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of a minimum of 3 minutes.
- 4) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
- 5) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

Headquarters
Oticon A/S
Kongebakken 9
DK-2765 Smørum
Denmark



SBO Hearing A/S
Kongebakken 9
DK-2765 Smørum
Denmark

244223US / 2022.09.09 / v2