Beltone **Solus Max**™

6 reasons why

Beltone Solus Max will improve your fitting experience



Download now on gnhearing.com

Hearing care makes a big difference to the lives of many. Provide your very best hearing care and spend time where it matters to you, with a versatile and easy fitting experience.

A smoother and faster experience

It's easier than ever to find what you need for your fittings, with 74% less clutter[1] on the home screen, and a new user-friendly discovery screen featuring intuitive navigation.

Do what you need in fewer clicks

With our smarter software, you can save more than a minute of time per fitting,[2] find settings with a 33-second faster click-through speed[1]and 62% fewer clicks[1] – and update both hearing aids simultaneously in 30% less time due to streamlined firmware updates.[3]

Complete verification testing

Get ready to experience our most comprehensive setup for standard verification and directionality testing. Complete all verification measurements while staying connected to the hearing aids, in 50% less time than before.[4]

Showcase the sound experience

The extended range of demo styles* lets your clients experience our technologies and benefits first-hand. Switch between technology levels with just two clicks.[5]

See all the options

Use the catalogue feature to quickly find and compare all Beltone hearing aids and features across technology levels.

Elevated onboarding experience

Quickly and easily access client report information directly in Noah, a client report session on the data logging screen, and personalised reports for your clients in their preferred language.

"I like it so much more. It's less verbiage, it's less overwhelming"

Hearing care professional, 1 year of experience, USA

"It's eye-catching and premium. The Connect button is now at the center. I can find it more easily"


Hearing care professional, 5 years' experience, China

1 Improvements compared to Beltone Solus Max 1.17. Technical data on file
2 Time saving of 1 minute and 19 seconds per fitting on average, compared to Beltone Solus Max 1.18. GN Proprietary data on file
3 Overall time saving of 30% compared to Beltone Solus Max 2.1. GN Proprietary data on file (2025)
4 50% time saving compared to Beltone Solus Max 2.1 (from 5.36 mins to 2.54 mins)
5 GN Proprietary data on file
* Flexible demo devices are available as microRIE, miniRIE, RIE312, RIE13 styles and Super Power BTE


Beltone™

6 easy steps to fit Beltone hearing aids


1. Launch Beltone Solus Max fitting software and put the devices in discovery mode by removing them from the charger or closing the battery door for non-rechargeable models. Click **Connect to Beltone Solus Max**.

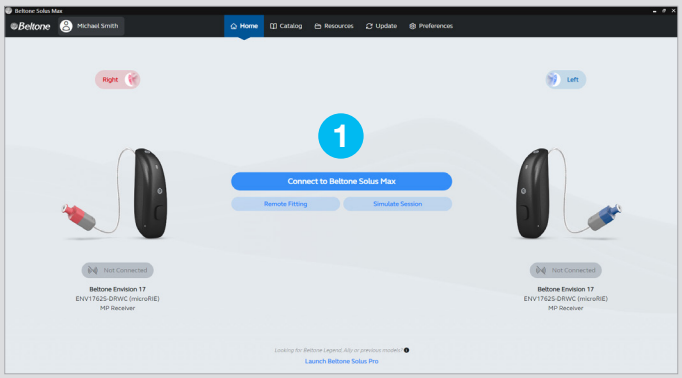


2. The hearing aids will appear in the **Select Instruments** section. Use the **Push to Assign** feature to assign to the respective right and left ear placement. Click **Continue** in the lower right of the screen. Select appropriate physical properties and click **Continue** to proceed through the connection process.

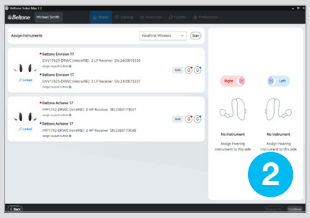


3. Calibrate the hearing aids on your client's ears to activate Feedback Eraser 2. Click **Continue** in the lower right of the screen. **Note**: Devices will automatically unmute once calibrated.

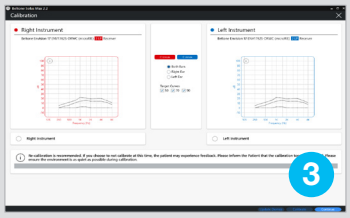

4. Navigate to **Patient > Profile** screen, customise the Prescribed First Fit selections as appropriate:
   a. Select Patient **Experience Level** based on previous experience with amplification and to maximise acceptance of hearing aids.
   b. Click **Fit Patient** in the lower right, or navigate to **Gain Adjustments**.

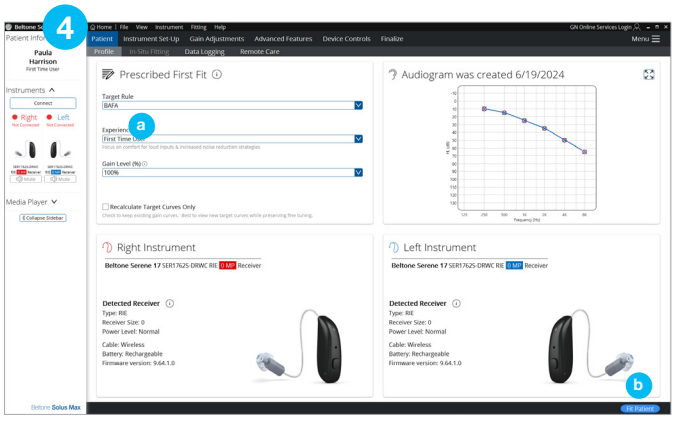

5. Fine-tune instruments using **Gain Adjustments**:
   a. Add or remove programs from this screen.
   b. Once programming is set, navigate to **Device Controls** to customise the manual push button, demonstrate beeps or pair accessories.

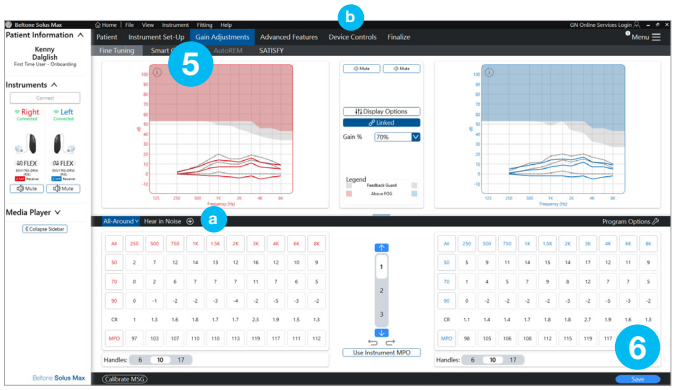

6. Click **Save** at the bottom right of the screen. Devices will reboot and play the start-up melody upon exiting the software.

Worldwide Headquarters
GN Hearing A/S
Lautrupbjerg 7
DK-2750 Ballerup
Denmark
Tel.: + 45 45 75 11 11
gnhearing.com

M201709GB-25.06-Rev.A

© 2025 GN Hearing A/S. All rights reserved. Beltone is a trademark of GN Hearing A/S.