Beltone **Solus Max**™

6 reasons why

Beltone Solus Max will improve your fitting experience



Download now on gnhearing.com

Hearing care is not just a job. It's something that makes a big difference in the lives of many. And as a hearing care professional your expertise is key to providing the best possible care. So, we want to provide brand-new and versatile software, so you can spend time where it matters to you.

Find what you want in fewer clicks

With a 33-second improvement to click-through speed[1], and a 62% click reduction[1] when navigating the home screen and preferences, we're making it easier for you to find what you need. We're bringing Beltone Solus Max into the modern world with a modern look. We've simplified our user interface, so there's now 74% less clutter[1].

Feature and technology comparison

It's now possible to compare features across different hearing aid technology levels. You can get an easy overview of what the different levels of technology can offer your clients.

Show – don't just tell

Introducing our new range of RIE (Receiver-in-Ear) demo styles that lets you change technology levels easily with 2 clicks[2], so you can let your clients hear and experience the difference advanced sound features make.

Personalised onboarding for your clients

Set up and print personalised onboarding reports for your clients in Beltone Solus Max. That way, they can see what programs you've given them, what the push buttons do, and so much more.

Easy access to counselling and updates

We've moved the preferences menu to the front, so you can access all of our updates to the software and firmware, as well as our many great counselling tools, such as how-to videos, straight from the home screen.

Fittings done your way

You can set up your favourite settings to become part of your standard fitting flow. Beltone Solus Max will also keep track of your fittings, and when you use the same settings again and again, it will recognise it, and propose it to become a standard part of your fitting flow.

"I like it so much more. It's less verbiage, it's less overwhelming"

Hearing Aid Fitter, 1 years experience, USA

"It's eye-catching and premium. The Connect button is now at the center. I can find it more easily"


Hearing Aid Fitter, 5 years experience, China

Beltone™


1 Beltone Solus Max 2.0 compared to version 1.17, technical data on file
2 GN Proprietary data on file

6 easy steps to fit Beltone hearing aids


1. Launch Beltone Solus Max fitting software and put the devices in discovery mode by removing them from the charger or closing the battery door for non-rechargeable models. Click **Connect to Beltone Solus Max**.

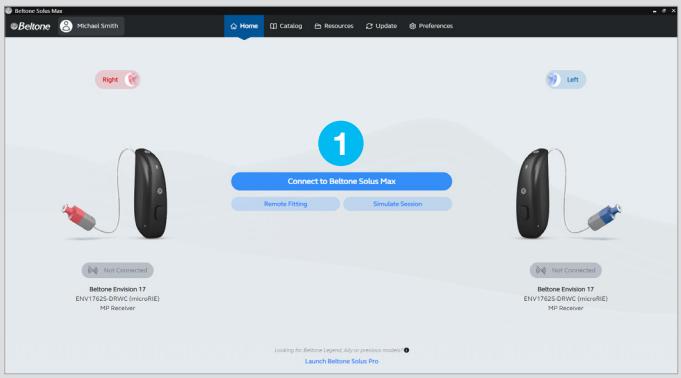


2. The hearing aids will appear in the **Select Instruments** section. Use the **Push to Assign** feature to assign to the respective right and left ear placement. Click **Continue** in the lower right of the screen. Select appropriate physical properties and click **Continue** to proceed through connection process.

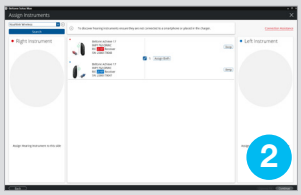


3. Calibrate the hearing aids on the patient's ears to activate Feedback Eraser 2. Click **Continue** in the lower right of the screen. **Note**: Devices will automatically unmute once calibrated.

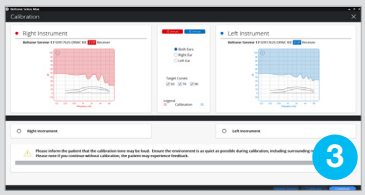

4. Navigate to **Patient > Profile** screen, customise the Prescribed First Fit selections as appropriate:
   a. Select Patient **Experience Level** based on previous experience with amplification and to maximise acceptance of hearing aids.
   b. Click **Fit Patient** in the lower right, or navigate to **Gain Adjustments**.

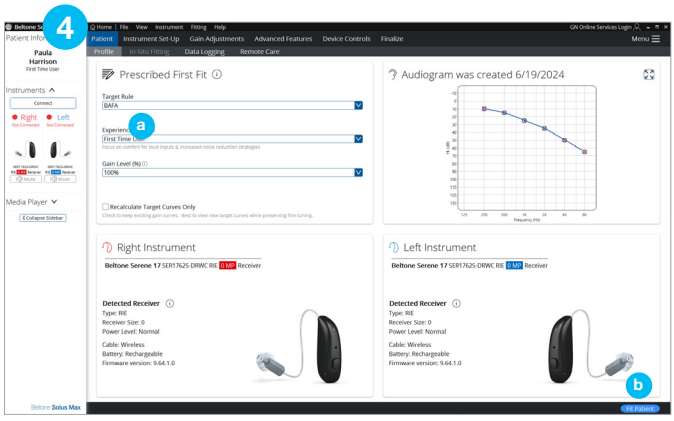

5. Fine-tune instruments using **Gain Adjustments**:
   a. Add or remove programs from this screen.
   b. Once programming is set, navigate to **Device Controls** to customise the manual push button, demonstrate beeps or pair accessories.

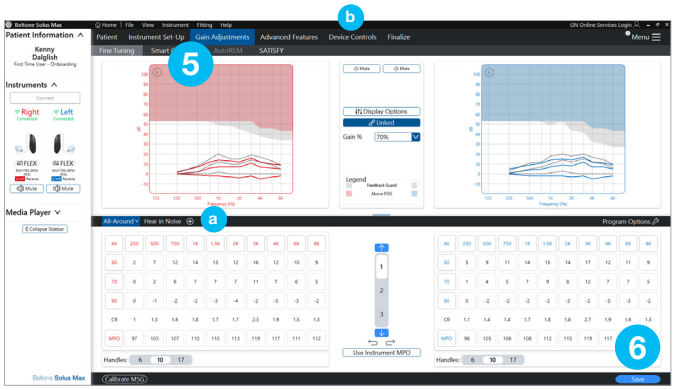

6. Click **Save** at bottom right of screen. Devices will reboot and play the start-up melody upon exiting the software.

Worldwide Headquarters
GN Hearing A/S
Lautrupbjerg 7
DK-2750 Ballerup
Denmark
Tel.: + 45 45 75 11 11
gnhearing.com

M102986GB-25.01-Rev.A

© 2025 GN Hearing A/S. All rights reserved. Beltone is a trademark of GN Hearing A/S.