
For clients who stream with an Android™ smartphone

Beltone wireless hearing aids enable direct streaming of audio from Android smartphones and tablets – the result of a unique collaboration with Google. Beltone was the first to offer this capability with direct streaming from Android devices to Beltone Amaze™ and now Beltone Boost Ultra™ hearing aids. Check if your client's Android device is compatible and how to get started streaming.

For information on which Android devices are compatible, visit: beltone.com/compatibility

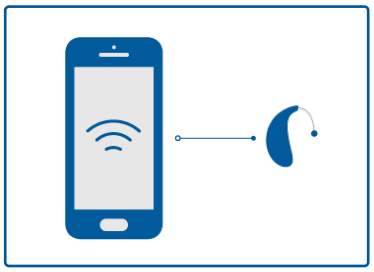

Our app can help you check Android streaming compatibility

1. Go to Google Play. In the Search field on top, enter "Beltone HearMax" and download the Beltone HearMax™ app.
2. Once downloaded, open the app and tap Get Started.
3. Follow the instructions on the screen to connect. You are now ready to go.

Compatibility scenarios

1 GOOD TO GO!
- Smartphone is compatible.
- Hearing aids are paired in Bluetooth® settings as Hearing Aid Devices

No red badges will be visible on:
- My Beltone menu
- Direct Audio Streaming menu point

2 MUST PAIR CORRECTLY
- Smartphone is compatible
- Hearing aids are not paired in Bluetooth settings as Hearing Aid Devices


Red badges will be visible on:
- Status and My Beltone menus
- Direct Audio Streaming menu point

Tap the menu point for guidance to pair correctly.

3 NOT COMPATIBLE
- Smartphone is not compatible

Red badges will be visible on:
- My Beltone menu
- Direct Audio Streaming menu point

Tap the menu point for information about compatibility requirements.

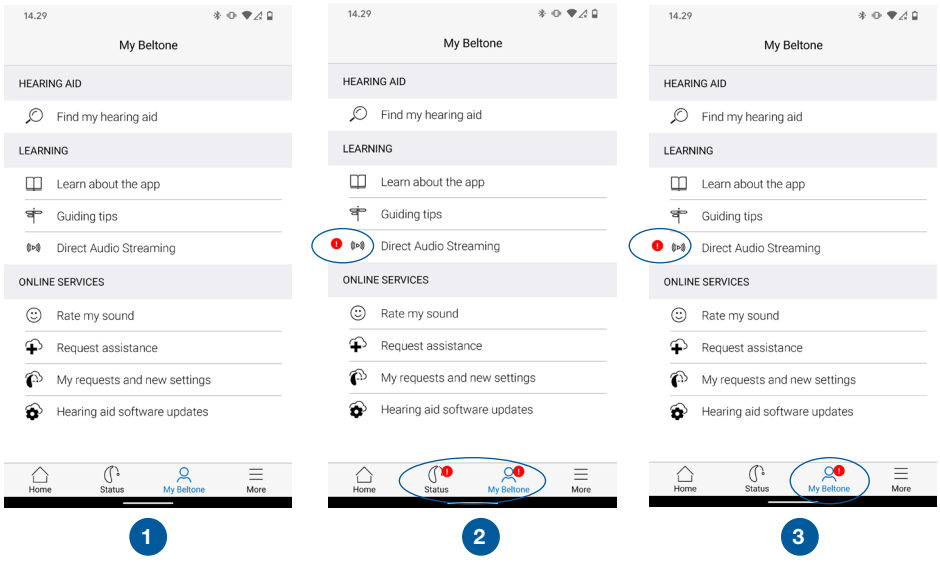

NOTE: The direct audio streaming requires minimum Android 10 and Bluetooth 5.0, and that Android Streaming to Hearing Aids has been enabled by the phone manufacturer. If direct Android streaming is not supported by your hearing aids or Android smartphone, you can use the Beltone Phone Link 2 to stream sound. The Beltone Phone Link 2 enables you to connect your Beltone hearing aids to any Bluetooth v. 2.0 Android device.

M102778GB-20.02-Rev.A

Beltone™

© 2020 GN Hearing Care Corporation. All rights reserved. Beltone is a trademark of GN Hearing Care Corporation.Android is a trademark of Google LLC. The Bluetooth word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.