Charger overview

ReSound offers rechargeable solutions with inductive charging and long battery life to keep your patients on the go.

How to match rechargeable hearing aids with charger

Each charger type is engineered to hold the hearing aid in the optimal position for charging reliability and efficiency. The charging form factor (CFF) relates to a specific ReSound hearing aid collection or style.

The CFF is indicated with a number (see chart below). The number is shown inside the charging bay and on the charger packaging.



	Premium charger	Standard charger	Desktop charger	Custom charger*
microRIE ReSound Sensia™, ReSound Sensia™ Bimodal, ReSound Vivia™ & ReSound Nexia™	8	8	8	-
miniRIE ReSound Savi™	9	-	9	-
BTE ReSound Nexia™	-	-	70	-
Power/High Power BTE ReSound Nexia™	-	-	80	-
Super Power BTE ReSound Enzo™ IA & ReSound Enzo™ IA Bimodal	-	-	90	-
ITC/ITE ReSound Sensia™ & ReSound Nexia™	-	-	-	C

*Custom charger includes custom-made inserts. Remind patients to bring their custom-made charger along with their custom hearing aids to all follow-up appointments.

Charger features


	ReSound Sensia, ReSound Viva & ReSound Nexia				ReSound Enzo IA	ReSound Savi	
	Premium charger	**Standard charger**	**Desktop charger**	**Custom charger**	**Desktop charger**	**Premium charger**	**Desktop charger**
	microRIE		microRIE/BTE/ Power BTE	ITE/ITC	Super Power BTE	miniRIE	
Inductive charging	•	•	•	•	•	•	•
Portable with onboard power bank	•					•	
Portable without power bank		•		•			
USB-C connector and cable	•	•	•	•	•	•	•
LED on the back	Shows the power bank battery status	Shows when the charger is connected to a power supply	Shows when the charger is connected to a power supply	Shows when the charger is connected to a power supply	Shows when the charger is connected to a power supply	Shows the power bank battery status	Shows when the charger is connected to a power supply
LED on the front	Shows battery level of the hearing aids*			Shows the hearing aids charging		Shows the battery level of the hearing aids*	


*Required to open the lid to locate LED

How to tell if the hearing aids are charging

Place the hearing aids in the charger and let the light indicator on the hearing aids guide you. For rechargeable Customs, the LED on the front of the charger will be your guidance.

You know that charging has started when the hearing aid lights blink slowly.

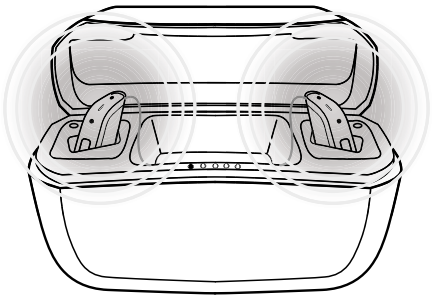

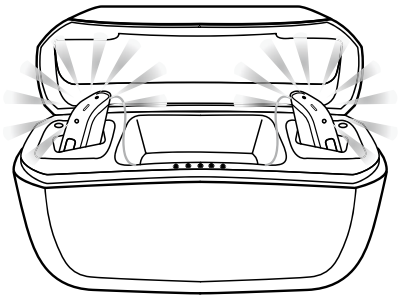

When fully charged, the hearing aid lights remain constant, not blinking.

The lights turn off when the hearing aids are removed from the charger.

Premium charger illustrated. The process is the same for the Standard and Desktop chargers.
For more information, please see the individual user guide for each charger.

Manufacturer according to FDA:

GN ReSound North America
5005 Dean Lakes Blvd.
Shakopee, MN 55379
USA
1-800-248-4327
pro.resound.com

© 2026 GN Hearing A/S. All rights reserved. ReSound is a trademark of GN Hearing A/S.

MK605546 Rev J 2026.06