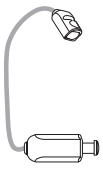
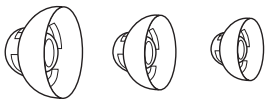
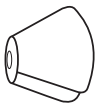
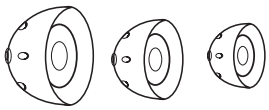
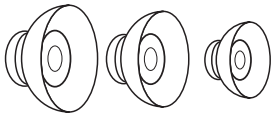
ReSound SureFit 2C RIE Receiver and Dome Overview

Multiple receiver and dome options enable ReSound Receiver-In-The-Ear (RIE) hearing aids to accommodate a wide range of hearing losses.

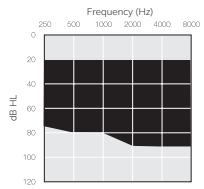
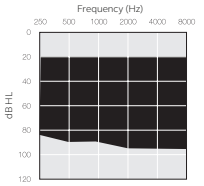
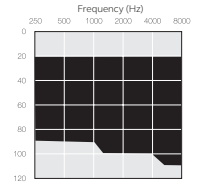
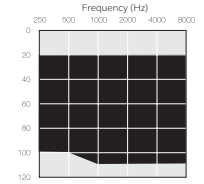
Receiver domes ReSound SureFit 2C receivers are compatible with SureFit 3 domes	LP: Low Power 	MP: Medium Power 	HP: High Power 	UP: Ultra Power 
OPEN (SM, MDM, LGE) Provides the patient with the least occluding experience. 	✓	✓	✓	UP receivers require encased molds
TULIP More occluding than open domes. Known to be more comfortable than other domes for some patients. One size option makes it a good choice for demos. 	✓	✓	✓	
CLOSED (SM, MDM, LGE) Occlusion similar to Tulip Domes but less occluding than Power Domes. 	✓	✓	✓	
POWER (SM, MDM, LGE) Most occluding dome option for the least risk of feedback. 	✓	✓	✓	
GN Wax Filters (Small, White)	✓	✓	✓	✓

ReSound SureFit 2C Receiver Guide

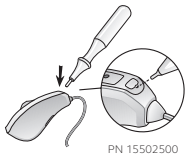
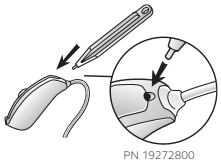
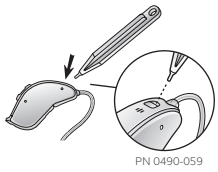
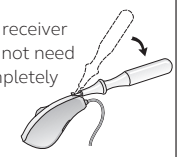
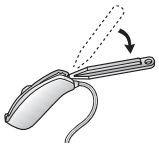
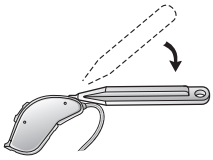
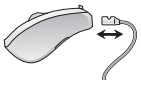
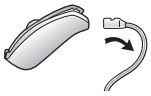
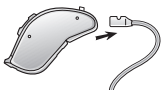
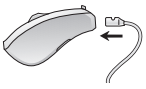
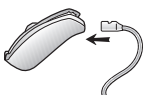
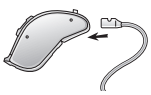
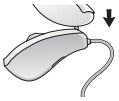
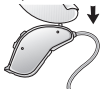


61 Lithium-ion | 61 312 | 62 13

Receiver Power Levels

	ReSound LiNX Quattro ReSound Key All 61 and 62 RIEs
LP **Low power receiver** 	
MP **Medium power receiver** 	
HP **High power receiver** 	
UP **Ultra power receiver** 	

Changing a SureFit 2C Receiver

ReSound LiNX Quattro ReSound Key* Rechargeable 61	ReSound LiNX Quattro ReSound Key 61	ReSound LiNX Quattro ReSound Key 62
Insert one end of the white tool into the opening of the clear receiver lock just above the top microphone.  PN 15502500	Insert the tip of the white tool into the small slot on the housing, behind the receiver.  PN 19272800	Insert the tip of the blue tool into the small slot in the receiver locking mechanism.  PN 0490-059
Using the white tool as a lever, lift the receiver lock gently until it protrudes from the hearing instrument. Note: the receiver lock does not need to be completely removed. 	Using the tool as a lever, swing the tool out and away from the hearing instrument to unlock the microphone cover. 	Lower the tool at the angle shown, and pull the receiver locking mechanism up. 
Pull the receiver out to remove it from the hearing instrument. 	Pull the receiver down and out to remove it from the hearing instrument. 	Pull the old receiver tube straight out to remove it from the hearing instrument. 
To insert a new receiver, push the new receiver tube until it is fully seated. 	To insert a new receiver, push the new receiver tube in until it is fully seated. 	To insert a new receiver, push the new receiver tube in until it is fully seated. 
Lock the receiver tube in place by pushing the clear receiver lock into the hearing instrument until it is flush.  Attach a dome or custom Micromold to the receiver.	Lock the receiver tube in place by clicking the locking mechanism closed.  Attach a dome or custom Micromold to the receiver.	Lock the receiver tube in place by firmly pushing down the locking mechanism until snugly in place.  Attach a dome or custom Micromold to the receiver.

*ReSound Key Rechargeable hearing instrument available in US only.

Manufacturer according to FDA:

GN ReSound North America
8001 E Bloomington Freeway
Bloomington, MN 55420
USA
1-800-248-4327
pro.resound.com

ReSound Government Services
8001 E Bloomington Freeway
Bloomington, MN 55420
USA
1-800-392-9932
gs.resound.com

Manufacturer according to Health Canada:

ReSound Canada
2 East Beaver Creek Road, Building 3
Richmond Hill, ON L4B 2N3
Canada
1-888-737-6863
pro.resound.com

© 2021 GN Hearing A/S. All rights reserved. ReSound is a trademark of GN Hearing A/S.

MK605403 Rev A 2021.01