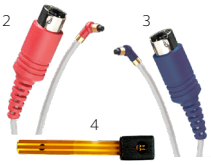
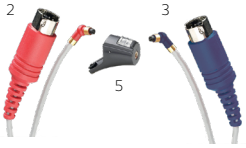
ReSound Smart Fit cable guide

Programming interface overview

	ITEM NUMBER	DESCRIPTION
1	23053200	Noahlink Wireless
2	21984000	CS44 – Red, Long
3	21983900	CS44 – Blue, Long
4	19533900	CS53 Flex Strip
5	19902500	PRG ADAPTER 2
6	19215000	Programming Kit, ENZO Q


	1  Noahlink Wireless	2 3 4  CS53 flex strip & CS44 socket cable (R&L)	2 3 5  Programming adapter & CS44 socket cable (R&L)	6 2  ENZO Q Programming adapter & CS44 socket cable (R&L)
ReSound Sensia™	RIE, Wireless Customs	CIC	-	-
ReSound Vivia™	RIE	-	-	-
ReSound Savi™	RIE, Wireless Custom	CIC	-	-
ReSound Nexia™	RIE, BTE, Wireless Customs	CIC	-	-
ReSound Enzo™ IA	SPBTE	-	-	-
ReSound ENZO Q™	BTE	-	BTE 88-DWHT	BTE 98-DWT
ReSound Key™	RIE, BTE, Customs	Customs, RIE 61-DRW, BTE 67-DW	RIE 62-DRW, BTE 77-DW, BTE 88-DWH	BTE 98-DW

NOTE: Strongly recommend to use fresh batteries for non-rechargeable models.

Cable connection for non-wireless programming

Custom models


With the battery door ajar, insert the flex strip into the slot beside the hinge pin. The black dot on the flex strip should be facing away from the battery door. The flex strip is fully inserted when the black dot abuts the hearing instrument case. Close the battery door once the flex strip is inserted.

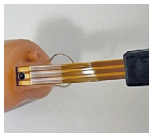



RIE 61-DRW

Locate the programming slot below the push-button. Open the battery door. With the black dot on the flex strip facing the push button on the hearing instrument, insert the flex strip into the slot and close the battery door.

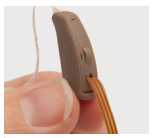



Mini BTE 67-DWT

Locate the programming slot below the microphone port. With the black dot on the flex strip facing the push button on the hearing instrument, insert the flex strip into the slot, until the black dot abuts the hearing instrument case.






RIE 62-DRW(T), BTE 77-DWT & BTE 88-DWHT

To connect the programming adaptor, hold it at a 45° angle and insert the small pins at the tip of the adaptor into the square cut-out on the underside of the hearing instrument. Once in place, move the adaptor in the direction of the battery compartment and then gently click it into place.



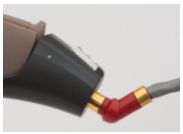

BTE 98 -DW(T)

[1] Locate and push out the battery door pin using the pin removal tool. **[2]** Pull out the battery door and **[3]** add the ENZO Q programming adapter (19215000) and secure it with the pin. **[4]** Insert a new size #675 battery and close the door. **[5]** Connect CS44 socket cable.

1. 
2. 
3. 
4. 
5. 

ReSound battery guide

	BEHIND-THE-EAR HEARING AIDS						
	RECEIVER-IN-EAR (RIE)			BEHIND-THE-EAR (BTE)			
BATTERY	Integrated Lithium-ion	#312	#13	Integrated Lithium-ion	#312	#13	#675
ReSound Sensia	SNx60S-DRWC	SNx61 -DRW	SNx62-DRW	-	-	-	-
ReSound Vivia	VIx60S-DRWC	VIx61-DRW	VIx62-DRW	-	-	-	-
ReSound Savi	SAx60-DRWC	SAx61-DRW	SAx62-DRW	-	-	-	-
ReSound Nexia	NXx60S-DRWC CX160S-DRWC	NXx61-DRW	NXx62-DRW	NXx77-DWC NXx88-DWC CX177-DWC	-	-	-
ReSound Enzo IA	-	-	-	EIx98-DWHC	-	-	-
ReSound ENZO Q	-	-	-	-	-	EQx88-DWHT	EQx98-DWT
ReSound Key	KEx61-DRWC*	KEx61-DRW	KEx62-DRW	-	KEx67-DW	KEx77-DW KEx88-DWH	KEx98-DW

*Available in US only

	CUSTOM HEARING AIDS				
BATTERY	Integrated Lithium-ion	#10A	#312	#312	#13
ReSound Sensia	SNxITE-DWC SNxITC-DWC	SNxCIC	SNxITC-DW	-	-
ReSound Savi	-	SAxCIC	SAxITC-DW	-	-
ReSound Nexia	NXxITC-DWC NXxITE-DWC	NXxCIC	NXxITC-DW	-	-
ReSound Key	-	KExCIC	KExITC-DW	KExITE-DW	KExITE-DW

INSTALLING NEW BATTERIES

To ensure successful programming, always use a fresh battery. Remove the protective foil from the battery and wait for one to two minutes before inserting the battery into the hearing aid to allow the battery to fully activate.

Manufacturer according to FDA:

GN ReSound North America

5005 Dean Lakes Blvd.
Shakopee, MN 55379
USA
1-800-248-4327
pro.resound.com

Manufacturer according to Health Canada:

ReSound Canada

2 East Beaver Creek Road, Building 3
Richmond Hill, ON L4B 2N3
Canada
1-888-737-6863
pro.resound.com

© 2026 GN Hearing A/S. All rights reserved. ReSound is a trademark of GN Hearing A/S.

MK604817 Rev AB 2026.06