
6 reasons why

ReSound Smart Fit™ will enrich your fitting experience

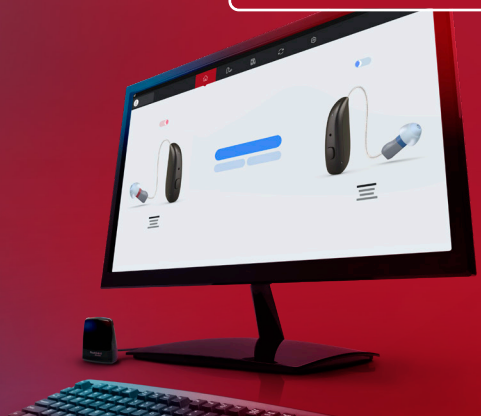

Download now on gnhearing.com

As a hearing care professional, you spend a lot of your day in fitting software.
You deserve the best experience. Enjoy smarter and more intuitive software, every time you fit.

Enjoy a more intuitive experience

The path to what you need has never been cleaner, simpler or faster. Experience 74% less clutter[1] on the home screen and explore a new revamped, user-friendly discovery screen.

Spend less time clicking

Save more than a minute per fitting,[2] with 62% fewer clicks[1] to find preferences and settings; 33-second faster[1] click-through speed; and 30% less time needed for firmware updates.[3]

Complete verification testing

We're introducing our first directionality presets for test box testing via our fitting software to verify directionality. Complete all verification measurements while staying connected to the hearing aids – in 50% less time than before.[4]

Show – don't just tell

Let your clients experience our advanced sound features and benefits first-hand with our great range of demo styles.* Switch between technology levels with just two clicks.[5]

Compare across technology levels

Quickly look up, check and compare all ReSound hearing aid models and features across technology levels in our catalogue feature.

Greater onboarding experience

Access client report information directly in Noah, a report session on the data logging screen, and personalised client reports in the preferred language.

We've prioritised the fundamentals

ReSound Smart Fit keeps on evolving and improving to meet the needs of hearing care professionals around the world, continuously enhancing your fitting software experience.

"I think it is visually very appealing. It's larger, it's easier to read, it's less cumbersome with all the stuff popping up."

Hearing Care Professional, 7 years' experience, USA

"Sometimes you need a quick reference on what the serial number is. It's nice that's available here."


Hearing Care Professional, 1 year of experience, USA

1 Improvements compared to ReSound Smart Fit 1.17. Technical data on file
2 Time saving of 1 minute and 19 seconds per fitting on average, compared to ReSound Smart Fit 1.18. GN Proprietary data on file
3 Overall time saving of 30% compared to ReSound Smart Fit 2.1. GN Proprietary data on file (2025)
4 50% time saving compared to ReSound Smart Fit 2.1 (from 5.36 mins to 2.54 mins)
5 GN Proprietary data on file


* Flexible demo devices are available as microRIE, miniRIE, RIE312, RIE13 styles and Super Power BTE

6 easy steps to fit ReSound hearing aids


1. Launch ReSound Smart Fit fitting software and put the devices in discovery mode by removing them from the charger or closing the battery door for non-rechargeable models. Click **Connect to ReSound Smart Fit**.

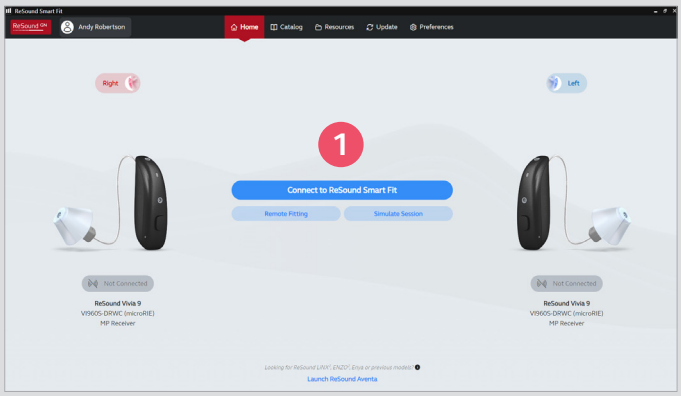


2. The hearing aids will appear in the **Select Instruments** section. Use the **Push to Assign** feature to assign to the respective right and left ear placement. Click **Continue** in the lower right of the screen. Select appropriate physical properties and click **Continue** to proceed through the connection process.

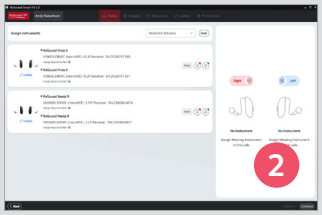


3. Calibrate the hearing aids on your client's ears to activate DFS Ultra III. Click **Continue** in the lower right of the screen. **Note**: Devices will automatically unmute once calibrated.

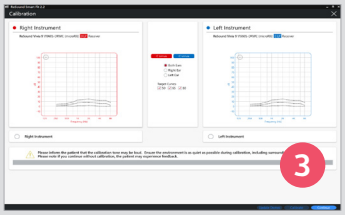

4. Navigate to **Patient > Profile** screen, customise the Prescribed First Fit selections as appropriate:
   a. Select Patient **Experience Level** based on previous experience with amplification and to maximise acceptance of hearing aids.
   b. Click **Fit Patient** in the lower right, or navigate to **Gain Adjustments**.

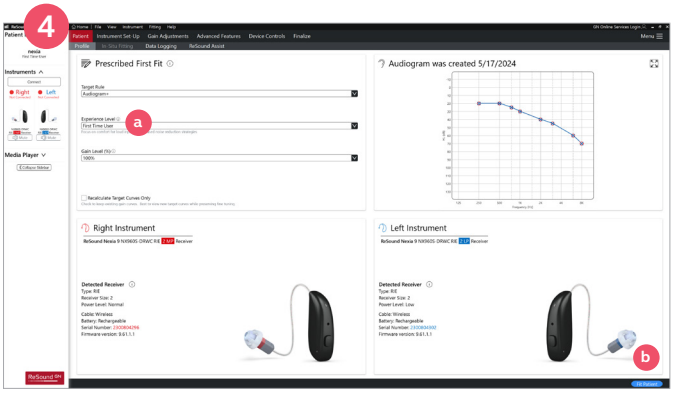

5. Fine-tune instruments using **Gain Adjustments**:
   a. Add or remove programs from this screen.
   b. Once programming is set, navigate to **Device Controls** to customise the manual push button, demonstrate beeps or pair accessories.

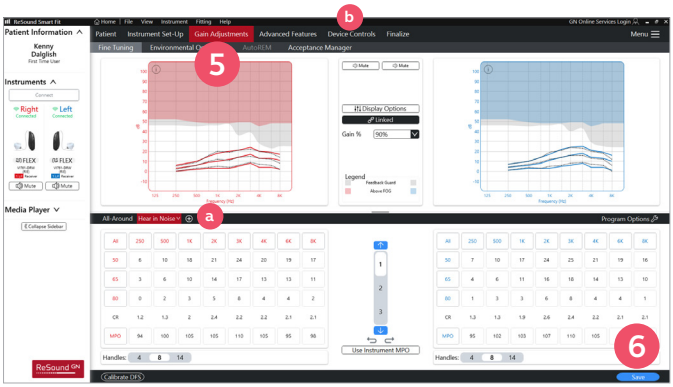

6. Click **Save** at the bottom right of the screen. Devices will reboot and play the start-up melody upon exiting the software.

Worldwide Headquarters
GN Hearing A/S
Lautrupbjerg 7
DK-2750 Ballerup
Denmark
Tel.: +45 4575 1111
gnhearing.com

CVR no. 55082715

M2017708GB-25.06-Rev.A

© 2025 GN Hearing A/S. All rights reserved. ReSound is a trademark of GN Hearing A/S.