Beltone Rely
Fitting Flow

A guide for professionals



Beltone **Rely**™

Beltone™



First fitting with Beltone Solus Max™ 1.10

This fitting guide gives an overview of how to fit Beltone Rely wireless hearing instruments with Beltone Solus Max 1.10. A Beltone Rely 463 DRWC wireless binaural fitting is depicted. Not all features described in this fitting guide are applicable for Beltone Rely fittings.

Beltone Solus Max 1.10 fitting software supports all Beltone Rely, Beltone Imagine™, Beltone Amaze™, Beltone Trust™ and Beltone Boost Max™ hearing instruments.

Prior to connecting to Beltone Solus Max

Rechargeable hearing instruments can only be connected through wireless interface and cannot be connected with cables. Ensure that the instruments are charged prior to connection.
Note: They cannot be connected while in the charger.


Insert Noahlink Wireless in a USB slot on the fitting PC and wait for the message that new hardware has been found. Then open Beltone Solus Max.

Launch Beltone Solus Max. For new fittings, Beltone Solus Max begins in the Launcher.

Launcher

Both Beltone Solus Pro 1 and Beltone Solus Max 1.10 are contained within one software installation. The Launcher will detect both Beltone Rely and legacy instruments, and then retrieve the software appropriate for the connected device.


1. Reboot the hearing instrument by opening and closing the battery door.
2. Click ‘Connect’ in the upper right of the Launcher screen.

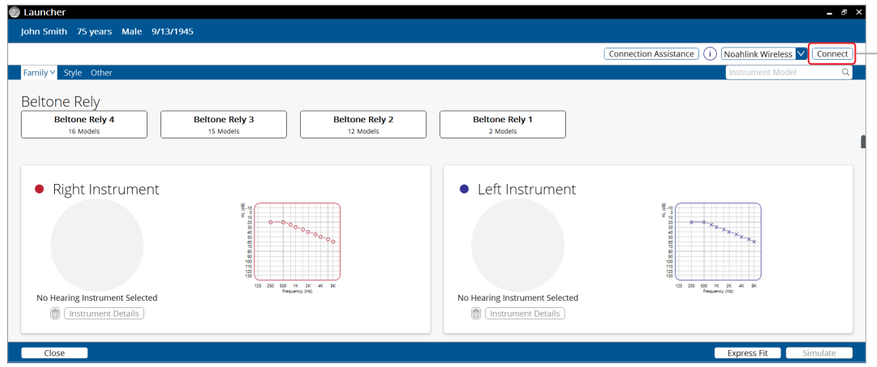

Select instruments

The hearing instruments will appear in the ‘Select instruments’ section.


1. When the hearing instruments to be programmed have been identified select them by checking the left box.
2. Assign the hearing instruments to the right and left ears by choosing the correct radio button for the right and left sides. When pressing “Beep Instrument”, the audio signal will be played. It is recommended to beep the instrument for confirmation that the instrument selected in the software is in fact the instrument being used in the fitting and on the correct side.
3. Continue in the lower left of the screen once the hearing instruments are selected and assigned.

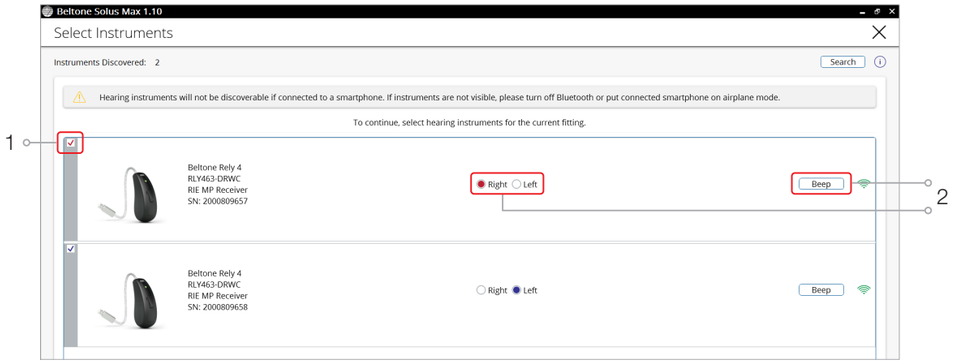


Connection flow

Once the hearing instruments have been selected the connection process will begin.


1. Select receiver power of the instrument.

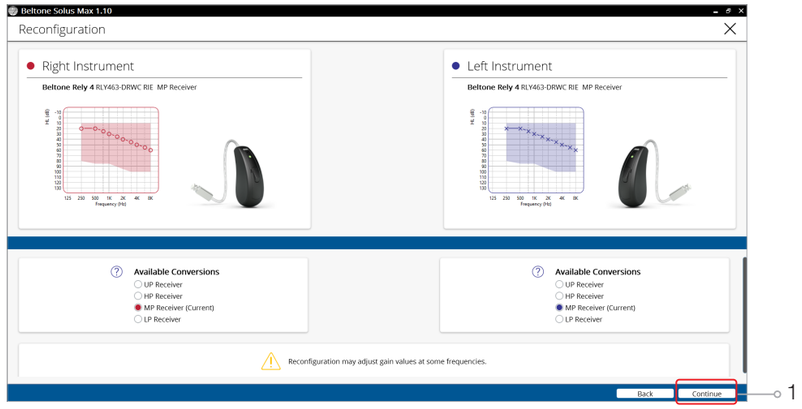


2. Select physical properties of the instrument.

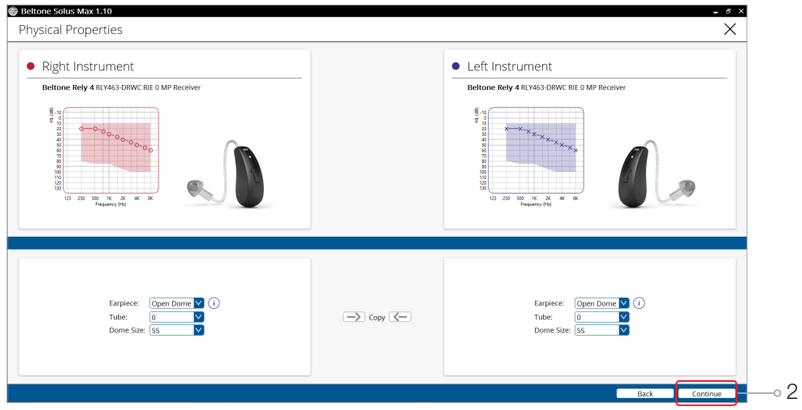

3. Calibrate instruments to activate Feedback Eraser.

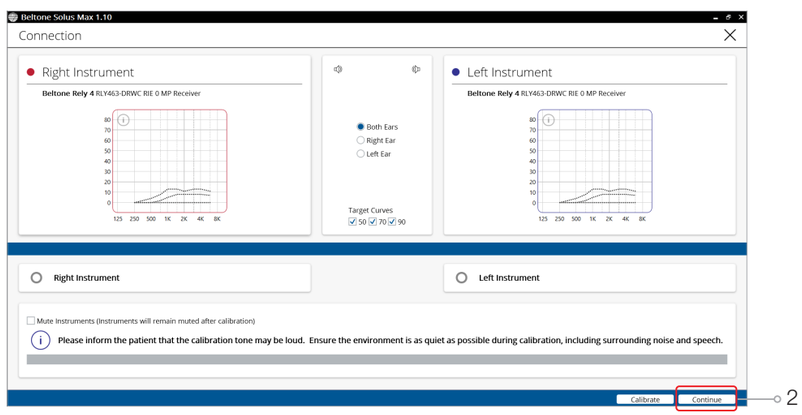


4. Verify connected hearing instruments. Click ‘Continue’ to proceed to the Patient Profile screen or pair accessories by clicking ‘Accessories’.

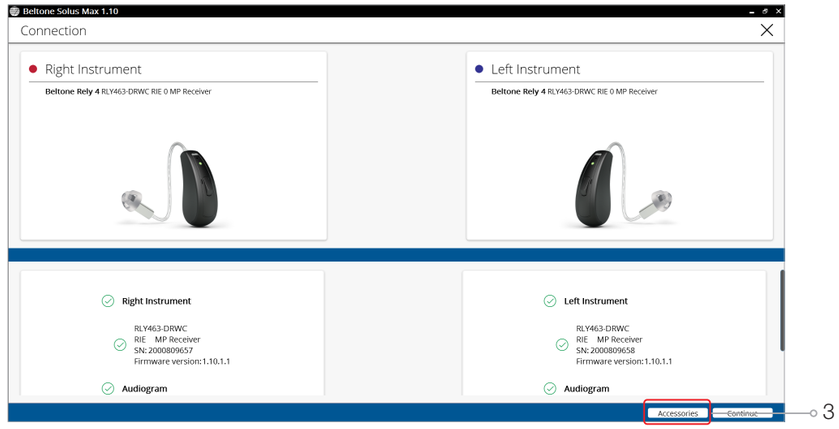

Patient profile screen

Located in the upper left corner of Patient Profile screen, Prescribed First Fit ensures an immediate personalisation of a new set of hearing aids. Based on a user’s audiogram, the selected target rule, gain level (%) and experience level; the Beltone Solus Max software will adjust the settings for Sound Cleaner Pro, Impulse Noise Reduction and Silencer to better fit the individual user.

1. Select the Patient Experience Level based on the amplification history.
2. Set Gain Level(%) for initial fit.
3. Click ‘Fit Patient’ button in the lower right corner to move to Fitting Screen.

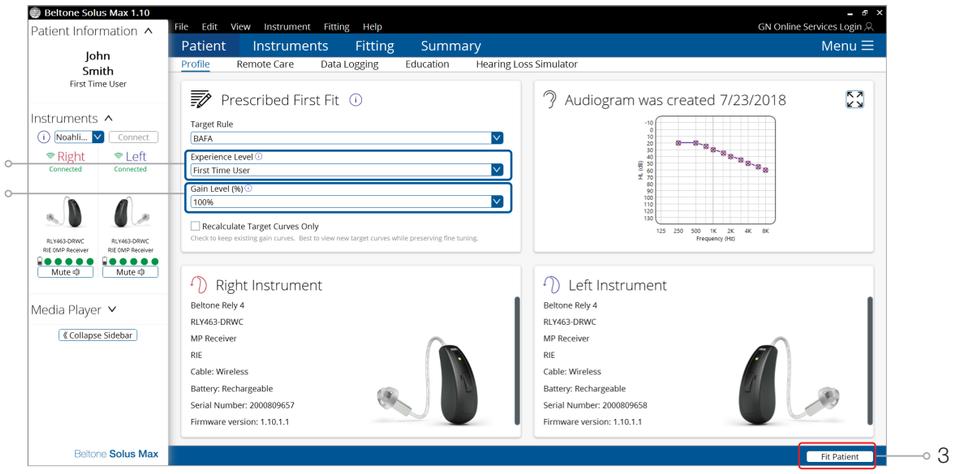


Note the Beltone Remote Care feature tab in the lower level navigation under the blue upper navigation. The new remote services feature can be found and activated in the “Remote Care” tab on the Patient screen.
If the patient being fitted is a candidate for Beltone Remote Care, the feature must be activated during the fitting in order to be active in the hearing instrument and appear in the Beltone Hear Max app.

Pairing wireless accessories

The Accessories screen can be accessed via the connection flow at the final step, at the bottom of the fitting screen, or in the Accessories tab of the Instruments navigation item.


1. Click ‘Pair’ beneath the accessory you wish to pair.

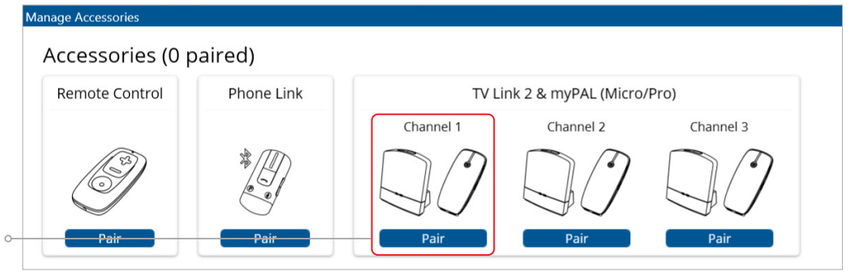

2. Once the searching pop-up appears, press the pairing button on the accessory with a small pointed objects such as a pen tip or paper clip. Follow the pairing instruction displayed.



3. When the accessory is successfully paired, the accessory image will appear in color with a remove button below.

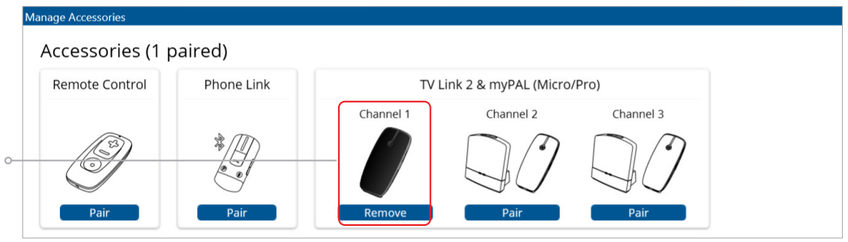

Fitting screen

On the Fitting screen, gain adjustments can be made as well as modification to advanced features.

1. Gain adjustments can be to one or several frequencies by selecting frequencies between 250 Hz and 6,000 Hz on the gain grid and then by clicking the up and down arrows found between the right and left ear grids. Gain can be adjusted for soft, average, and loud input levels which alters average the compression ratios noted above the MPO.


2. Access advanced features and to instrument controls is also found on the Fitting screen. The tabs below the blue navigation bar will provide access these tools.

3. Additional tools and features can be found in the upper right menu.

4. Programs are listed in the program tabs. To select a different program select arrow in tab to view program options. Select the ‘+’ to add additional programs. Program tools can be found in the tab on the far right.

5. Click the ‘Save’ button to save fitting to hearing instruments and database.

6. Wireless accessories can be paired both outside and inside the software. To pair inside the software select ‘Pair Accessories’ to access the Accessories screen. Be sure that the accessory is charged or plugged in, and turned on during the pairing process.

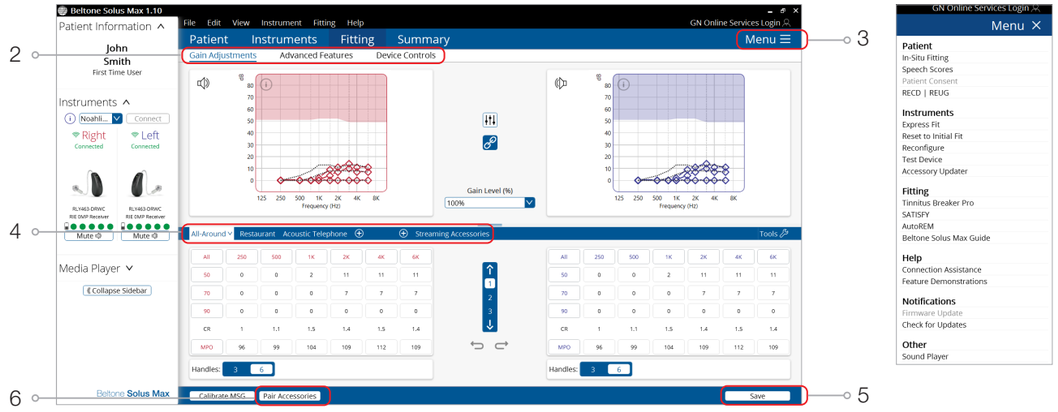


Note: For some users, fine tuning may be needed to meet certain gain targets.

Advanced features

Advanced features are prescribed per program and can be modified from their default settings selecting from the drop down box or by moving the slider to the desired value. The defaults are set based on the program environment which is selected in the program tab.

Directionality

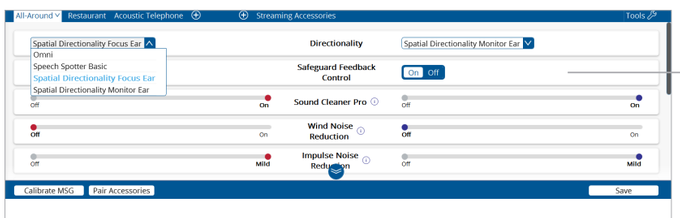

Note the information icons for a short explanation of the feature. Some information text may be clickable and direct to an in depth explanation and video.


Spatial Directionality: The Focus (directional) ear and Monitor (omnidirectional) ear are prescribed based on the audiogram and any speech audiometry data that may be available.

Fixed Hypercardoid: Fixed hypercardoid directionality is the simplest directional option and provides a non-varying hypercardioid pattern. In this setting, signals coming from behind and the sides are reduced in order to enable the user to better concentrate on signal from the front. The directional characteristics are constant and static. If selected, this type of directionality is always ‘on.’

Speech Spotter Basic: Speech Spotter Basic Directionality is similar to Fixed Directionality, but Adaptive Directionality dynamically changes the directional pattern to reduce the loudest sounds from the rear. This is preferable to Fixed Hypercardoid Directionality for users often in environments with multiple, moving, or simultaneous noise sources.

Omni: Omnidirectionality means that the microphones are equally sensitive to sounds from all directions. Omni is often preferred for sound quality when listening in quiet environments.

Smart Beam: Smart Beam Directionality is an advanced form of Adaptive Directionality that allows the width of the directional beam to be set in the fitting software, in effect making the area of focus larger or smaller.

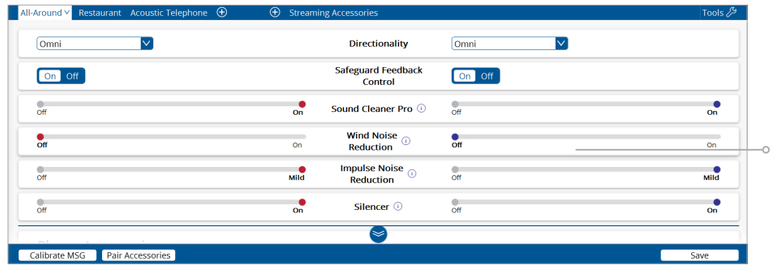

By scrolling down or pulling up the lower part of the Advanced Features screen all the features available for the connected hearing aid can we viewed on the screen.

Feedback Eraser with WhistleStop

This is the feedback control. It is activated when the devices have been calibrated for Feedback Eraser. The default setting is Mild. Adjust the setting to a stronger level (Moderate or Strong) if concerns for feedback arise. The Music Mode setting is a Mild setting, typically used for music programs or if the patient experiences feedback only when tonal sounds occur (for example, from a musical instrument).

Wind Noise Reduction

Wind Noise Reduction applies noise reduction specifically for windy situations. The default setting is Off.
If a patient is experiencing too much wind noise. If a patient is experiencing too much wind noise, turn on.

Impulse Noise Reduction


Sudden impulse sounds can be quite uncomfortable for those using hearing instruments.
Impulse Noise Reduction detects the loud impulse sound and attenuate it instantaneously.

Safeguard Feedback Control

Safeguard Feedback Control is a feedback control setting that is activated automatically before the hearing instrument has been calibrated for Sound Cleaner Pro. Once Sound Cleaner Pro has been calibrated, this setting is hidden, and cannot be reselected unless the hearing instrument has been restored to factory settings. Safeguard Feedback Control is designed for use when a patient is unable to return to the clinic for the fitting, and the hearing instrument is instead shipped to the patient from the fitter. It can be turned off if the fitter does not believe the patient will need any feedback control, in the case of very low gain settings. It applies a mild level of feedback cancellation, but as it has not been calibrated for the patient's ear, it may not be as effective at controlling feedback as a Mild setting of Feedback Eraser with WhistleStop.

Phone Accessories

Features related to the streamed sound from a Phone Link, or a smart device that is capable of direct streaming, are located at the lower portion of each environmental program tab. Streamed sound will inherit the gain settings of the environmental program that is being used during the time of streaming. While the phone accessory features can be adjusted within the individual program tabs, it is important to note that settings are applied across programs, therefore phone accessory features changes made in one program will also be reflected in others. Please note when viewing Data Logging the time spent streaming from the phone will be collected within the manual program that is being used during streaming.

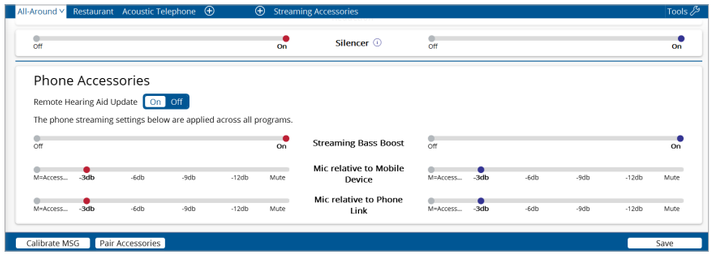

Completing a follow-up fit with Beltone Solus Max 1.10


Prior to connecting to Beltone Solus Max

Note: Rechargeable Beltone Rely hearing instruments can only be connected through wireless interface and cannot be connected with cables. Ensure that the instruments are charged prior to connection. Note that they cannot be connected while in the charger.

Insert Noahlink Wireless in a USB slot on the fitting PC and wait for the message that new hardware has been found. Then open Beltone Solus Max.

It is not possible for Noahlink Wireless to connect to hearing instruments that are actively connected to a smart device via Bluetooth. Ensure that the hearing instruments do not have an active Bluetooth connection to a smart device. This can be done by having the user turn off their smart device or disable the Bluetooth during the visit. If an active connection to Bluetooth is detected, Beltone Solus Max will provide a reminder to deactivate the connection with the smart device.

Launch Beltone Solus Max. For follow-up fittings, the software that was used in the previous fitting will launch.

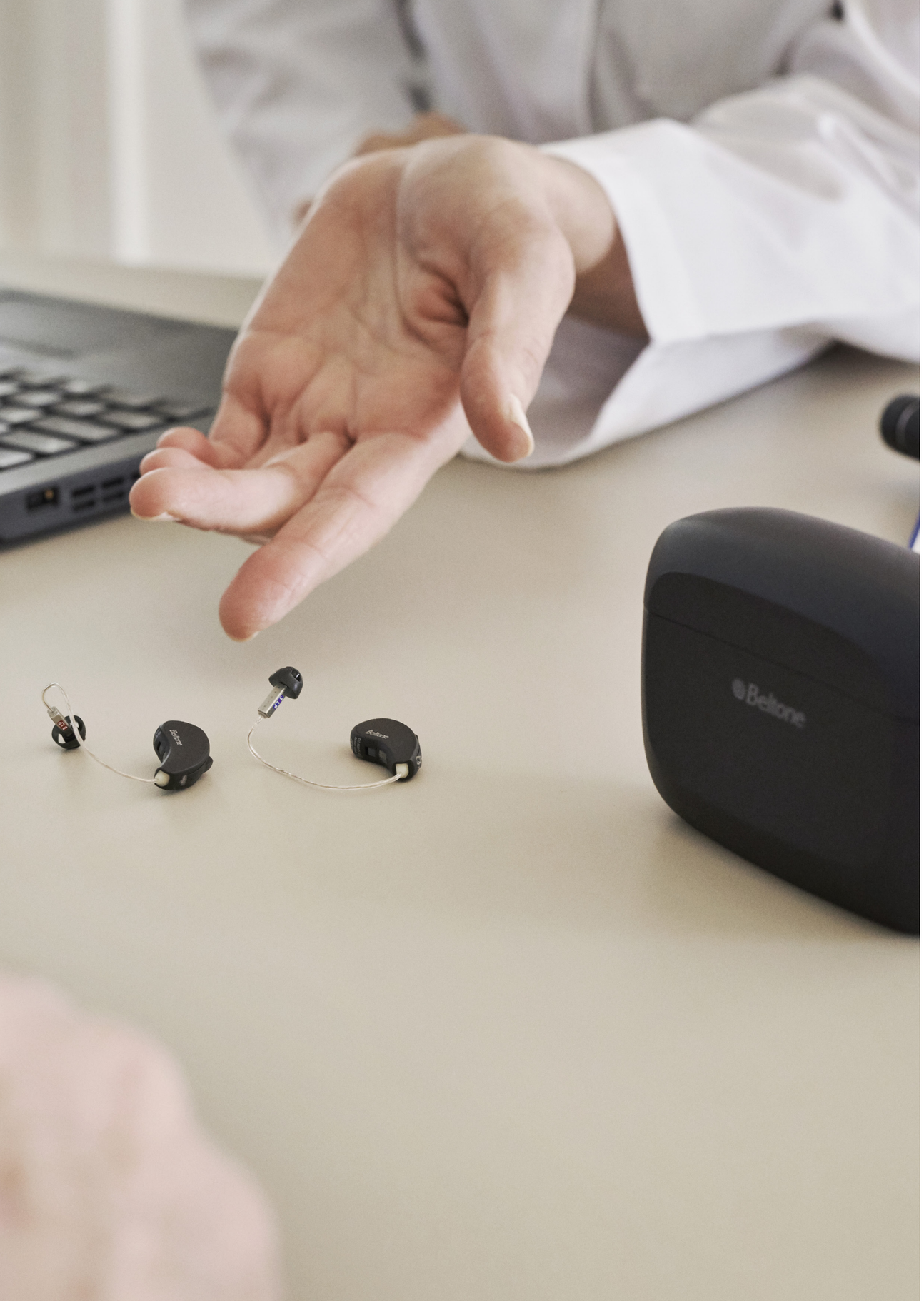

Upgrading a patient from legacy hearing instruments to Beltone Rely


Beltone Solus Pro 1

If the patient has an existing hearing instrument fit in Beltone Solus Pro, then Beltone Solus Pro 1 will open when you begin a new fitting.


1. Click the ‘Go to Solus Max’ button at the bottom left of the screen to close Beltone Solus Pro and arrive at the Solus Max Launcher.

2. Follow steps in the ‘First Fitting’ section of this guide

Beltone Rely hearing instruments cannot be fitted in Beltone Solus Pro software.

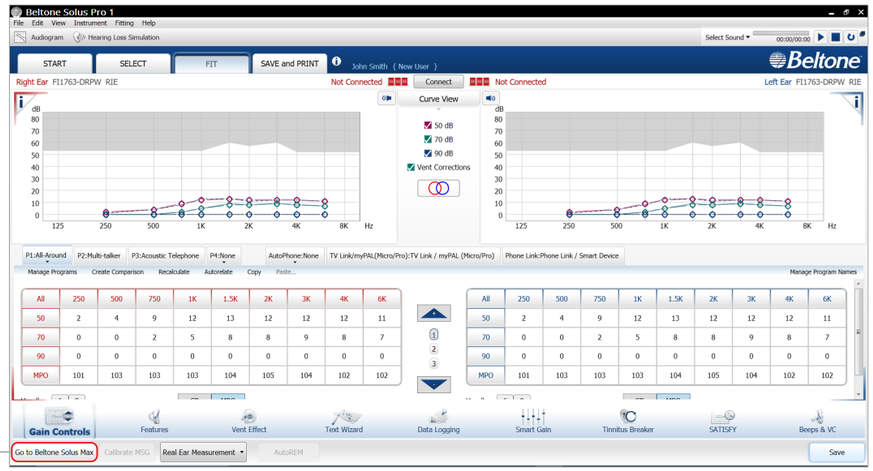

1

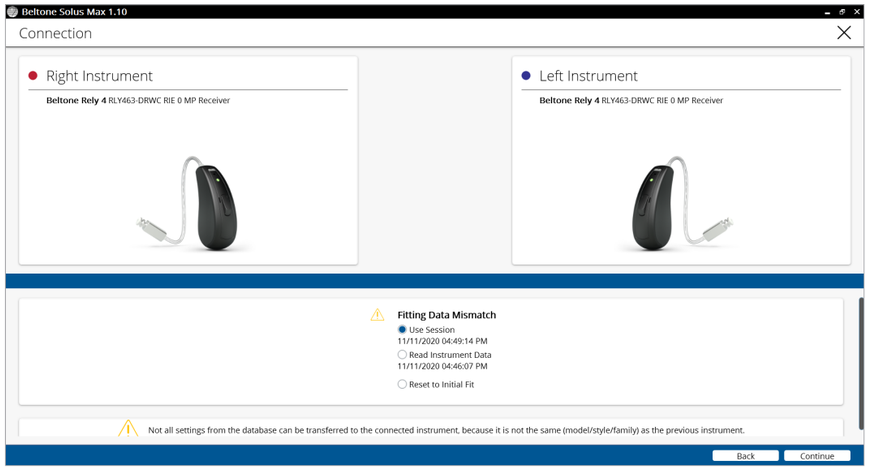


Note: To transfer setting from legacy fitting to Beltone Amaze fitting, select ‘use session’ in the fitting data mismatch screen of the connection flow. To start a fresh fitting choose ‘reset to initial fit’.

Replacing one hearing instrument of a pair

Connecting during a follow-up

When returning for a follow-up visit, the software will open to the last saved session. Depending on the saved preferences, the session will open to either the Patient Profile screen or the Fitting Gain Adjustments screen.

1. Open and close the battery door to reboot the hearing instrument.
2. Select Noahlink Wireless and click ‘Connect’ in the collapsible side bar panel.
3. The hearing instruments will appear in the ‘Select Instruments’ section. When the hearing instrument to be programmed has been identified, select it by checking the left box.
4. The hearing instrument that has already been programmed for this patient will already be assigned to the right or left side. When the other instrument is selected, it will automatically be assigned to the opposite side. It is recommended to beep the instrument for confirmation that the instrument selected in the software is in fact the instrument being used in the fitting and on the correct side.
5. Continue in the lower left of the screen once the instruments are selected and assigned.
6. Continue to connection flow and follow steps in the ‘First Fitting’ section of this guide.

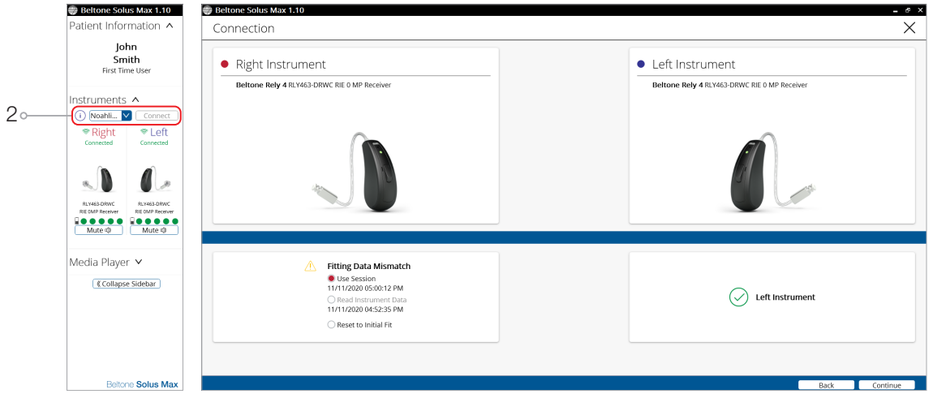

Note: To transfer setting from session to replaced instrument, select ‘use session’ in the fitting data mismatch screen of the connection flow. To start a fresh fitting choose ‘reset to initial fit’.

Beltone. Helping the world hear better

Beltone was founded on the act of helping a friend to enjoy life more. Since 1940, we have provided knowledge, tools, service and training to the professionals servicing the hearing impaired. Beltone makes quality hearing care available to more people. To do so, we develop technically optimal hearing solutions without ever losing focus on the individual's needs. It is our belief that no individual should be denied the possibility of an improved life with better hearing.

Worldwide Headquarters
Beltone A/S
Lautrupbjerg 7
DK-2750 Ballerup
Denmark
Tel.: + 45 45 75 11 11
Web: www.beltone.com
www.hearing.tv

Beltone™

© 2020 GN Hearing Care Corporation. All rights reserved. Beltone is a trademark of GN Hearing Care Corporation. Apple, the Apple logo and iPhone are trademarks of Apple Inc, registered in the U.S. and other countries. Android is a trademark of Google Inc.

M201512GB-20.11-Rev.A